


# PyMOL User's Guide

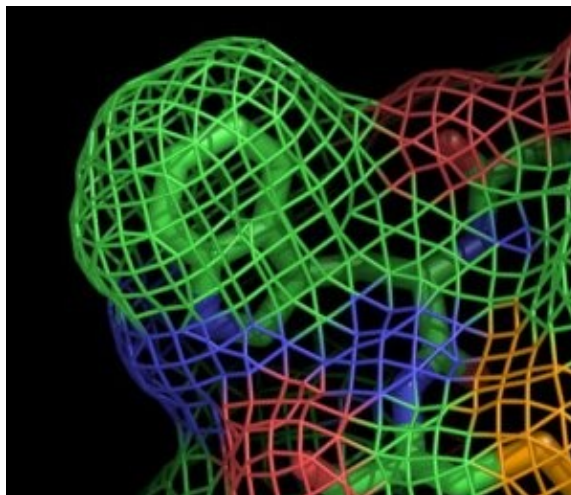

written by  
**Warren L. DeLano, Ph.D.**  
with assistance from  
**Sarina Bromberg, Ph.D.**

Copyright © 2004  
DeLano Scientific LLC  
All Rights Reserved.

# Table of Contents

<b><u>Copyright Notice and Usage Terms</u></b> .....	<b>1</b>
<u>Copyright Notice</u> .....	1
<u>Terms of Usage for the PyMOL User's Manual</u> .....	1
<u>Trademarks</u> .....	1
<b><u>Preface</u></b> .....	<b>2</b>
<u>Why PyMOL?</u> .....	2
<u>Words of Caution</u> .....	2
<u>Strengths</u> .....	3
<u>Weaknesses</u> .....	3
<b><u>Introduction</u></b> .....	<b>4</b>
<u>Welcome to PyMOL!</u> .....	4
<u>Is PyMOL Free Software?</u> .....	4
<u>Yes, but...</u> .....	4
<u>The DeLano Scientific Mission</u> .....	4
<b><u>Installation</u></b> .....	<b>6</b>
<u>Windows</u> .....	6
<u>Recommendations</u> .....	6
<u>Minimal System Requirements</u> .....	6
<u>Python-Free Installation</u> .....	6
<u>Python-Dependent Installation</u> .....	6
<u>MacOS X</u> .....	6
<u>Recommendations</u> .....	7
<u>Minimal Requirements</u> .....	7
<u>If you use Fink</u> .....	7
<u>If you do not use Fink</u> .....	7
<u>Linux and Unix</u> .....	8
<u>System Requirements</u> .....	8
<u>Dependency-Free Approaches</u> .....	8
<u>Dependency-Based Approaches</u> .....	9
<b><u>Getting Started with Mouse Controls</u></b> .....	<b>10</b>
<u>Launching</u> .....	10
<u>Using the Mouse</u> .....	10
<u>Using a Command Line</u> .....	10
<u>PyMOL's Windows</u> .....	11
<u>The Viewer Window</u> .....	11
<u>The External GUI Window</u> .....	12
<u>Loading PDB Files</u> .....	13
<u>Manipulating the View</u> .....	13
<u>Basic Mouse Control</u> .....	13
<u>Virtual Trackball Rotation</u> .....	14
<u>Moving Clipping Planes</u> .....	16
<u>Changing the Origin of Rotation</u> .....	16
<u>Getting Comfortable</u> .....	17

# Table of Contents

<b><u>Getting Started with Commands</u></b>	<b>18</b>
<u>Recording Your Work (Optional)</u>	18
<u>Loading Data</u>	18
<u>Manipulating Objects</u>	19
<u>Atom Selections</u>	19
<u>Coloring Objects and Selections</u>	21
<u>Turning Objects and Selections On and Off</u>	22
<u>Changing Your Point of View</u>	23
<u>Saving Your Work</u>	23
<u>Scripts and Log Files</u>	24
<u>png Files</u>	24
<u>Session Files</u>	25
<u>Command-Line Shortcuts</u>	25
<u>Command Completion using TAB</u>	26
<u>Filename Completion using TAB</u>	26
<u>Automatic Inferences</u>	26
<u>Other Typed Commands and Help</u>	27
<b><u>Command Syntax and Atom Selections</u></b>	<b>28</b>
<u>Syntax</u>	28
<u>Selection-expressions</u>	28
<u>Named Atom Selections</u>	29
<u>Single-word Selectors</u>	31
<u>Property Selectors</u>	31
<u>Selection Algebra</u>	34
<u>Atom Selection Macros</u>	35
<u>Calling Python from within PyMOL</u>	37
<b><u>Cartoon Representations</u></b>	<b>38</b>
<u>Background</u>	38
<u>Accessibility</u>	38
<u>Pretty and Correct</u>	38
<u>Customization</u>	41
<u>Cartoon Types</u>	41
<u>Fancy Helices</u>	44
<u>Secondary Structure Assignment</u>	45
<b><u>Ray-Tracing</u></b>	<b>46</b>
<u>Important Settings</u>	46
<u>Saving Images</u>	47
<u>png</u>	47
<b><u>Stereo</u></b>	<b>48</b>
<u>Introduction</u>	48
<u>Supported Stereo Modes</u>	48
<u>Crosseye Stereo</u>	48
<u>Walleye Stereo</u>	48
<u>Hardware Stereo</u>	48

# Table of Contents

<b><u>Stereo</u></b>	
<u>Generating Stereo Figures</u>	48
<b><u>Movies</u></b>	<b>49</b>
<u>Concepts</u>	49
<u>States and Frames</u>	49
<u>Important Commands To Know</u>	49
<u>load</u>	49
<u>mset</u>	49
<u>mdo</u>	50
<u>mmatrix</u>	50
<u>Simple Examples</u>	50
<u>Complex Examples</u>	51
<u>Previewing Ray-traced Movie Images</u>	51
<u>cache frames</u>	51
<u>mclear</u>	51
<u>Saving movies</u>	51
<u>mpng</u>	52
<b><u>Advanced Mouse Controls</u></b>	<b>53</b>
<u>Picking Atoms and Bonds</u>	53
<u>Example Usage of the "pk" Atom Selections</u>	53
<u>The "lb" and "rb" Selections</u>	53
<u>Conformational Editing</u>	54
<b><u>Crystallography Applications</u></b>	<b>55</b>
<u>Crystal Symmetry</u>	55
<u>load</u>	55
<u>symexp</u>	55
<u>Electron Density Maps</u>	56
<u>load</u>	56
<u>isomesh and isodot</u>	56
<b><u>Compiled Graphics Objects (CGOs) and Molscript Ribbons</u></b>	<b>57</b>
<u>Introduction</u>	57
<u>Molscript Ribbons</u>	57
<u>load</u>	57
<u>Using Molscript</u>	57
<u>Creating Compiled Graphics Objects</u>	58
<u>CGO Reference</u>	59
<u>load cgo</u>	59
<b><u>Callback Objects and PyOpenGL</u></b>	<b>61</b>
<u>Introduction</u>	61
<u>Example</u>	61
<u>load callback</u>	62

# Copyright Notice and Usage Terms

## Copyright Notice

The PyMOL User's Manual is Copyright © 1998–2004 DeLano Scientific LLC, San Carlos, California, U.S.A. All Rights Reserved.

## Terms of Usage for the PyMOL User's Manual

**This manual is NOT free.** It is a *PyMOL Incentive Product* created to help you use the program while also generating recurring sponsorship for the project. This manual is made available for evaluation via the "honor" system: **You may evaluate this manual for a continuous period of up to one year without obligation.** If you wish to continue using this document beyond the end of the evaluation period, then you must become a sponsor of the project by purchasing a PyMOL license and a subscription to maintenance and support from DeLano Scientific LLC (<http://www.delanoscientific.com>).

Of course, if you are willing to sponsor the project today, then please don't wait a full year to start. The sooner your sponsorship comes in, the sooner we can apply it to improve the software and documentation!

Existing PyMOL subscribers may use this manual for no additional cost. However, subscribers who do not renew their subscription upon expiration must discontinue use of this and all other PyMOL Incentive Products. Though we have no direct means of enforcing this, we ask, in recognition of our declared scientific mission, that you honor the trust placed in you.

PyMOL users who are unable to sponsor the project by purchasing a PyMOL license and maintenance subscription are welcome to use Open–Source versions of PyMOL and any free documentation that can be found on the internet.

## Trademarks

PyMOL, DeLano Scientific, and the DeLano Scientific Logo are trademarks of DeLano Scientific LLC. Macintosh is a registered trademark of Apple Computer Inc., registered in the U.S. and other countries. Windows is a registered trademark of Microsoft Corporation in the U.S. and other countries. Linux is a trademark of Linus Torvalds. Unix is a trademark of The Open Group in the U.S. and other countries. MolScript is a trademark of Avatar Software AB. All other trademarks are the property of their respective owners.

*This chapter last updated June 2004 by Warren L. DeLano, Ph.D.*

Copyright © 2004 DeLano Scientific LLC. All rights reserved.

# Preface

## Why PyMOL?

PyMOL is one lone scientist's answer to the frustration he encountered with existing visualization and modeling software as a practicing computational scientist.

Anyone who has studied the remarkable complexity of a macromolecular structure will likely agree that visualization is essential to understanding structural biology. Nevertheless, most researchers who use visualization packages ultimately run up against limitations inherent in them which make it difficult or impossible to get exactly what you need. Such limitations in a closed-source commercial software package cannot be easily surmounted, and the same is still true for free programs which aren't available in source form.

Only open-source software allows you to surmount problems by directly changing and enhancing the way software operates, and it places virtually no restrictions on your power and opportunity to innovate. For these reasons, we believe that **open-source software is an intrinsically superior research product** and will provide greatest benefit to computer-assisted scientific research over the long term.

Launched over Christmas break in December 1999, PyMOL was originally designed to: (1) visualize multiple conformations of a single structure [trajectories or docked ligand ensembles] (2) interface with external programs, (3) provide professional strength graphics under both Windows and Unix, (4) prepare publication quality images, and (5) fit into a tight budget. All of these goals have since been realized. Although PyMOL is far from perfect and lacks such desirable features such as a general "undo" capacity, it now has many useful capabilities for the practicing research scientist. We hope that you will find PyMOL to be a valuable tool for your work, and we encourage you to let us know what ideas you have for making it even better.

## Words of Caution

**About the Manual:** This version of the manual has been updated for PyMOL version 0.86 (January 2003) but is still quite rough. Prepare yourself for omissions, errors, and potentially obsolete information. Make an informed decision to use the PyMOL manual at your own risk. Understand that this same caution applies to the program as a whole — you shouldn't be using PyMOL if you aren't willing to troubleshoot problems and take the initiative on the mailing list in order to discover solutions.

**About the program:** PyMOL was created in an efficient but highly pragmatic manner, with heavy emphasis on delivering powerful features to end users. Expediency has almost always taken precedence over elegance, and adherence to established software development practices is inconsistent. PyMOL is about getting the job done now, as fast as possible, by whatever means were available. PyMOL succeeds in meeting important needs today, but we view it as merely an initial step in a promising direction.

In time, we hope that we and others will follow by creating PyMOL-like software platforms which meet the needs of users but also provides the design rigor and code quality necessary to enable broad participation of outside developers. Though PyMOL will undoubtedly continue to expand and improve over the next decade, we expectd that its long term impact will primarily be to inspire other development efforts having more time and resources, and which will undoubtable achieve greater heights.

That isn't to say that you can't find good things about PyMOL's internal design. Indeed, we believe that there are many successful and instructive aspects to the program. However, we just hope to appropriately calibrate your expectations with respect to the code you will find if you with to "dive under to hood". Though the

program is Open-Source, it is best thought of as a dense, semi-opaque tool, best extended through Python rather than as a C coding environment in which to embed new technologies.

## Strengths

- **Cross-Platform.** A single code base supports both Unix, Macintosh, and Windows, using OpenGL and Python and a small set of Open-source external dependencies.
- **Command-Line and GUI Control** Real world applications require both.
- **Atom Selections.** Arbitrary logical expressions facilitate focused visualization and editing.
- **Molecular Splits/Joins.** Structures can be sliced, diced, and reassembled on the fly and written out to standard files (i.e. PDB).
- **Movies.** Creating movies is as simple as loading multiple PDB files and hitting play.
- **Surfaces.** As good if not better than Grasp, and mesh surfaces are supported too.
- **Cartoon Ribbons.** PyMOL's cartoons are almost as nice as Molscript but are much easier to create and render.
- **Scripting.** The best way to control PyMOL is through reusable scripts, which can be written in the command language or in Python.
- **Rendering.** A built-in ray tracer gives you shadows and depth on any scene. You also render externally.
- **Output.** PNG files output from PyMOL can be directly imported into PowerPoint.
- **Conformational Editing.** Click and drag interface allows you to edit conformations naturally. Sculpting allows the molecule to adapt to your changes.
- **Expandability.** The PyMOL Python API provides a solid way to extend and interface.

## Weaknesses

- **User Interface.** Development has been focused on capabilities, not on easy-of-use for new users.
- **Documentation.** Only recently has any documentation become available.
- **Object-Orientation.** There is a single monolithic, functional API.
- **Electrostatics.** PyMOL is not yet a replacement for Delphi/Grasp.
- **No Mechanics Engine** Although PyMOL sports potent molecular editing features, you can't yet perform any "clean-up".

*This chapter last updated January 2004 by Warren L. DeLano, Ph.D.*

# Introduction

## Welcome to PyMOL!

Over the years, PyMOL has become a capable molecular viewer with support for animations, high-quality rendering, crystallography, and other common molecular graphics activities. It has been adopted by many hundreds (perhaps even thousands) of scientists spread over thirty countries. However, PyMOL is still very much a work in progress, with development expected to continue for years to come.

## Is PyMOL Free Software?

### Yes, but...

PyMOL is Copyrighted software that is Free for all parties to use, modify, and redistribute. Because of PyMOL's unusual status, you can be confident that the time you invest today in learning the package will provide you with long term utility no matter where your career happens to take you. You will never be required to pay software license fees in order to use Open-Source PyMOL or to share it with others who might find it useful.

Nevertheless, PyMOL is not free to develop, document, maintain, and support. If you decide to adopt the package, then you are asked and expected to contribute to the project in some manner. Although such contributions may take a variety of forms, most PyMOL supporters choose to sponsor the project by purchasing (usually through their school or employer) a license and a renewable subscription to maintenance and support. All such contributions are entirely voluntary since we have intentionally abandoned the usual means of compelling compliance. Instead, we depend on your *free will* to provide vital funding for development and to cover other necessary expenses. In return, we provide specific incentives (called *Incentive Products*) as a reward for helping to fund the project. Example incentive products include this manual, extra features, enhanced platform-specific binary versions, and various other conveniences

Please take this request seriously. If you value PyMOL, then it is clearly in your interest to sponsor it. To find out how to donate or to purchase a license, visit the PyMOL web site at <http://www.pymol.org>

## The DeLano Scientific Mission

DeLano Scientific LLC is a private vision-centered software company which owns, develops, and supports the PyMOL package. Our mission as a commercial entity is to create highly effective tools for scientific research and to distribute them as broadly as possible while still succeeding as a healthy business. As a "boot-strapped" company, DeLano Scientific is not beholden to any outside investors who would insist upon maximum returns on investment. Thus, we have the rare privilege of being able to place Scientific and Medical Progress ahead of Profit in our hierarchy of values.

We have chosen a free and open-source approach for PyMOL because we believe this strategy will have the greatest positive impact on humanity. Visualization is key part of understanding the nature of life at the molecular level, and powerful visualization tools need to be universally available to all students and scientists if we are to make rapid progress in biomedical research. If PyMOL is successful, then we hope to expand the scope of our endeavors to meet other critical research needs in related areas.



Growth of the DeLano Scientific will depend entirely on the willingness of PyMOL users to adopt, nurture, and advocate for our volitional approach to software funding. Eventually, we hope to evolve into a major provider of scientific software for biomedical research and be distinguished by the quality, openness, and accessibility of our products, the trusting and nonexploitive relationships we form with our customers, and our willingness to work with all parties in advancing scientific software technologies.

*This chapter last updated June 2004 by Warren L. DeLano, Ph.D.*

# Installation

## Windows

### Recommendations

- Windows 2000 or XP.
- A late-model 3D OpenGL compatible graphics accelerator card from nVidia, ATI, 3Dlabs or similar.
- 512 MB RAM (768 MB or 1 GB preferred).
- 3 Ghz Pentium 4 processor or similar.

### Minimal System Requirements

- Windows 98 and ME, or later. PyMOL will not run on Windows 95 and NT.
- 3D OpenGL compatible graphics accelerator card.
- 256 MB RAM.
- 500 Mhz Pentium 3 processor.

Unless you have prior experience with Python, we recommend installing a version of PyMOL which does not require an external Python interpreter. Avoid versions of which contain "-py21", "-py22", "-py23" or similar in the filename.

### Python-Free Installation

1. Download the ".zip" format archive. For example,

`pymol-0_90-bin-win32.zip`

2. Extract the .zip file using WinZip (Windows XP can open .zip files directly).
3. Double click on the "Setup" or "Setup.exe" icon in the folder.
4. Answer the questions which follow.

You can now launch PyMOL from the Start menu.

### Python-Dependent Installation

If you already have Python installed and wish to use PyMOL with that interpreter, the process is virtually identical. The only difference is that you need to download a version of PyMOL which matches your desired Python version in the filename. For example:

`pymol-0_90-bin-win32-py22.zip`

`Python-2.2.2.exe`

would work with  
available from  
<http://www.python.org>.

## MacOS X

## Recommendations

- Mac OS 10.2.x or 10.3.x.
- Dual 2.0+ Ghz G5 system.
- GeForce4 or Radeon 9x00 OpenGL accelerator.
- 1 GB of RAM.

## Minimal Requirements

- Mac OS X 10.2.x
- Single 833 Mhz G4 system (will run on less, but performance is poor).
- 3D OpenGL graphics acceleration.
- 512 MB of RAM (1 GB recommended).

## If you use Fink

PyMOL is part of the Fink ports collection.

```
sudo -s  
apt-get pymol install
```

should be sufficient to get a functioning instance on your system. However, it may not be the most recent version. We also highly recommend installation of Apple's X Server, which enables PyMOL to access your accelerated graphics hardware.

## If you do not use Fink

### Option 1: MacPyMOL

At the request of various Macintosh users, as well as Apple itself, we have created **MacPyMOL**, a special native Aqua version of PyMOL for the Macintosh. The latest version of MacPyMOL can be downloaded from <http://delsci.com/macpymol>.

However, note that this version is an Incentive Product only available to PyMOL sponsors (but students and teachers are exempt from this requirement). For more information on MacPyMOL, contact [support@delanoscientific.com](mailto:support@delanoscientific.com).

### Option 2: PyMOL X11 Hybrid

This version of PyMOL includes a native Aqua-based OpenGL window and an X11-based Tcl/Tk external GUI (graphical user interface). Before launching the PyMOL X11 Hybrid, you must have Apple's X11 server installed and launched. The advantage of using this version is that it is fully compatible with Open-Source PyMOL, and does not require Fink. However, unlike MacPyMOL, this version does not support direct export of QuickTime movies. This binary build is free, but not Open-Source.

1. Download the "pymol-0\_XX-bin-osx-x11-hybrid.dmg.gz" compressed disk image.
2. Extract the archive and mount the disk image.
3. Copy the "PyMOLX11Hybrid" folder to the main Applications folder on your hard disk.

You can then launch PyMOLX11Hybrid by double-clicking on the PyMOLX11Hybrid icon.

# Linux and Unix

## System Requirements

- 3D OpenGL graphics acceleration.

There are a several different ways to install PyMOL on Linux. Please consult the [PyMOL Web Site](#) for additional information.

## Dependency-Free Approaches

These do not require installing any other packages in a privileged location on your system. All you need to do is download a "tar"-ball appropriate for your system, such as the following:

- pymol-0\_93-bin-linux-libc6-i386.tgz (for Linux)
- pymol-0\_93-bin-irix65-r10k.tgz (for SGI)
- pymol-0\_93-bin-solaris8-sun4u.tgz (for Solaris)

issue the following commands

```
gunzip < pymol-...-bin-...-.tgz | tar -xvf -  
./setup.sh
```

which will install the program, and then

```
./pymol.com
```

will launch PyMOL. You may then want to make a symbol link for this file to **~/bin/pymol** for easy launching.

```
ln -s $PWD/pymol.com ~/bin/pymol
```

***Install a minimal dependency binary build.***

***Compile PyMOL from source along with the "ext" dependencies distribution.***

Because the installation process is often subject to change, please see the INSTALL file from the current distribution for detailed instructions. In summary,

1. Download, extract, configure, and compile the external dependencies.
2. Download and extract the current PyMOL source distribution.
3. Create a symbolic link from the external dependencies to "ext" in the PyMOL directory.
4. Configure compilation by copying and modifying a "Rules.make" from the setup directory to reflect your system.
5. Run "make" to build pymol.
6. Create a pymol.com specific to your installation location.

You should be able to launch PyMOL by running pymol.com. I usually symbolic link this file into my "bin" directory as "pymol".

## Dependency–Based Approaches

You must install the following packages on your system

```
python (2.x), tcl (8.x), tk (8.x), libpng (1.x), zlib (1.x),  
glut (3.x), glut-devel (3.x), pmw*, and numeric* (numpy)
```

(\* = not required for RPM packages.)

You then have several choices:

### ***Using RedHat binary packages (RPMs).***

```
rpm -i pymol-0.90-1.rh73.py22.i386.rpm
```

### ***Using Python's distutils to compile and install PyMOL as a standard Python module.***

```
python setup.py build      (as a user)  
python setup.py install    (as root)  
python setup2.py install   (as root)
```

You can now run PyMOL with `./pymol.com`.

### ***Using Makefiles with preinstalled system dependencies.***

Because the installation process is often subject to change, please see the INSTALL file from the current distribution for detailed instructions. In summary,

1. Download and extract the current PyMOL source distribution.
2. Configure PyMOL by copying and modifying a "Rules.make" from the "setup" directory to reflect your system.
3. Run "make" to build pymol.
4. Create a pymol.com specific to your installation location.

You should be able to launch PyMOL by running `pymol.com`, and it may be convenient to add a symbolic link from this file into your "bin" directory as "pymol".

*This chapter last updated June 2004 by Warren L. DeLano, Ph.D.*

# Getting Started with Mouse Controls

## Launching

### Using the Mouse

#### *On Windows:*

Click on the **Start** menu, follow it to **Programs** (or **All Programs** on Windows XP), and then release the mouse on **PyMOL**.

#### *On Mac OSX (native version)*

Double-click on the PyMOL icon in the **Applications** folder on your main hard drive

### Using a Command Line

Various command line options can be included under both Windows and Unix to automatically open files and launch scripts. See "launching" in the reference manual for more information on these options.

#### *On Windows:*

At the command prompt, issue:

```
c:\program files\delano scientific\pymol\pymolwin.exe
```


If PyMOL is installed somewhere nonstandard, then use the correct drive letter and path.

#### *On Unix, Linux, and MacOS X (Fink version)*

If you installed using a package such as an RPM, then there is a good chance that "**pymol**" is already in your path. If not, then edit **pymol.com** in the PyMOL distribution and make sure **PYMOL\_PATH** points to the actual location of the distribution. Enter **./pymol.com** to start pymol. You will probably want to create a link "**pymol**" from this file in to a "bin" directory in your path so that you can launch the program anywhere by simply entering "**pymol**".

# PyMOL's Windows

PyMOL normally starts with two windows: The Viewer Window and the External (Tcl/Tk) GUI Window.

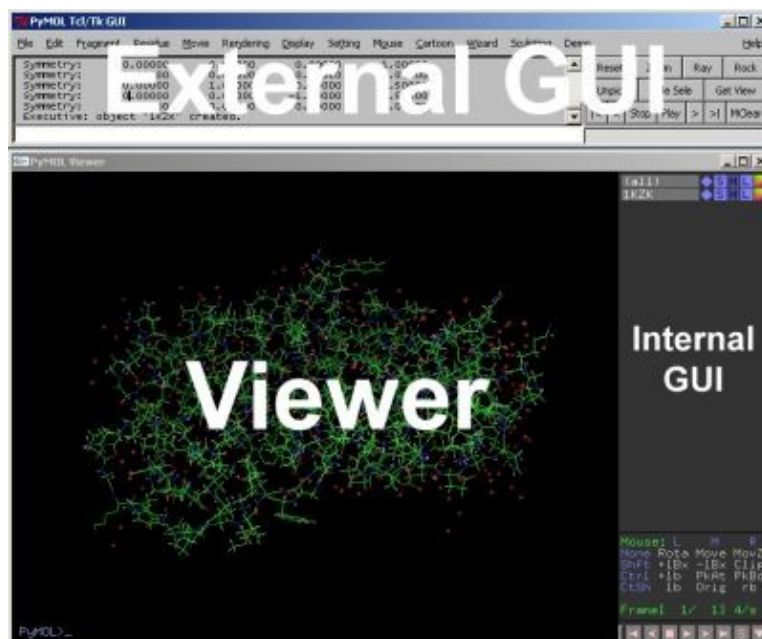


PyMOL's two windows.

GUI is an abbreviation for Graphical User Interface, which usually consists of menus, buttons, text boxes, and other familiar gadgets. By default, PyMOL actually has two GUI's: (1) an "Internal" GUI which appears inside the Viewer Window, and (2) an "External" GUI which appears inside of its own window. The reasons for this are boring and technical, but know that both GUI's will eventually be unified into a single interface in the future.

## The Viewer Window

The PyMOL Viewer represents the heart of the PyMOL system. This is a single OpenGL window where all 3D graphics are displayed and where all direct user interaction with 3D models takes place.

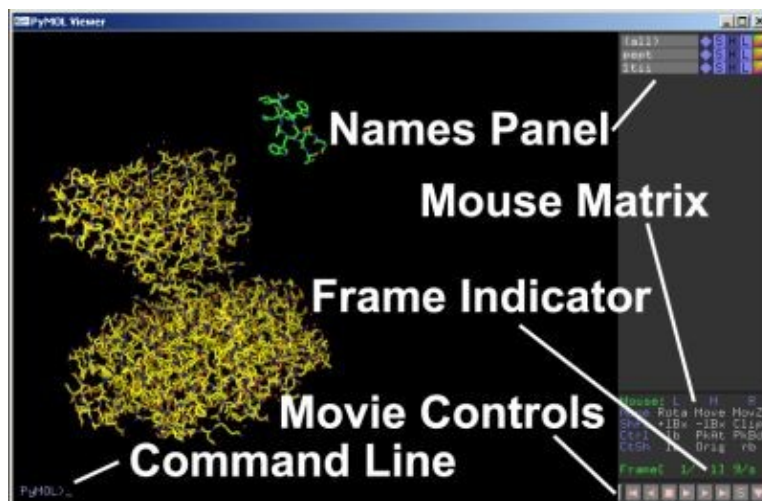


PyMOL Viewer window with Internal GUI enabled (Default).

The Internal GUI contained within this window (right) allows you to perform actions on specific objects and specific atom selections. From top to bottom, it contains an object list, a mouse button configuration matrix, a frame indicator, and a set of "VCR"-like controls for working with movies.

The Viewer also contains a command line (bottom) which can be used to enter PyMOL commands. It is also possible to view PyMOL text output in the Viewer window. you can **hit the ESC key anytime to toggle between text and graphics mode inside the Viewer window.**

The PyMOL Viewer can be run all by itself, and it provides the complete capability of the PyMOL core system. If desired, the Command line and Internal GUI can be disabled. Many tasks can be made easier and more efficient through use of standard menus and controls. For the most part, such gadgets are currently found in an External GUI window.

## The External GUI Window

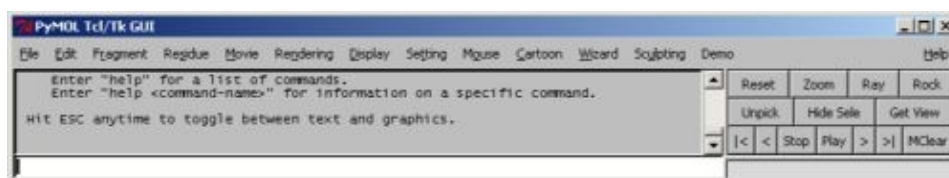

The default Tcl/Tk External GUI included with PyMOL.

By default, PyMOL comes with a single external GUI window which provides a standard menu bar, an output region, a command input field, and a series of buttons. One important advantage of the external GUI window is that **standard "cut and paste" functions for text will only work within the External GUI**, and not within in the PyMOL Viewer. Furthermore, **you must use Ctrl-X, Ctrl-C, and Ctrl-V to cut, copy, and paste** because a standard Edit menu has not yet been implemented.

**Notes For Developers:** External GUIs are the foundation for modularity and customizability in the PyMOL system. These windows constitute independent processes (or threads) which can control the behavior of PyMOL, and potentially interact with other programs. They are completely customizable at the Python scripting level, and mutiple external GUIs can exist at once (within the restrictions of Tkinter and wxPython).

External GUIs communicate with PyMOL through the Python API (Application Programming Interface). Those of you who want to link up you own programs with PyMOL should generally use a separate external GUI window to control the interaction, rather than changing internal PyMOL code. That way the programs will continue to work together even while development on each program proceeds independently. The internal GUI and all external GUI windows can be enabled and disabled using simple command line options (see reference for "launching").



## Loading PDB Files

### Using the External GUI Menu

The default external GUI provides a standard **Open...** item in the **File** menu which you can use to select the file you wish to open.


## Using Commands

## SYNTAX

```
load <filename>
```

### EXAMPLE

```
load test/dat/pept.pdb
```

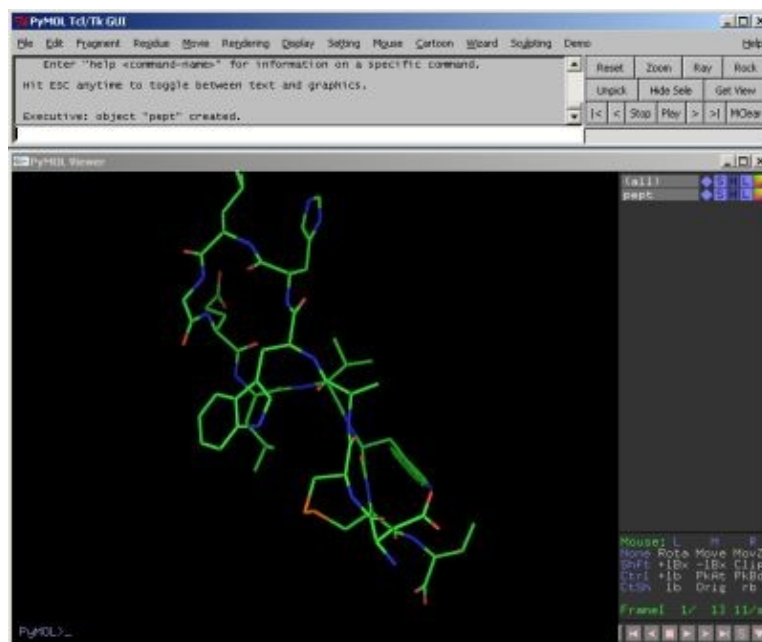


PyMOL after loading a PDB file.

## Manipulating the View

In PyMOL, the mouse is the primary control device, and keyboard modifier keys (SHIFT, CTRL, SHIFT+CTRL) are used in order to modulate button behavior. **A three button mouse is required for effective use of PyMOL**, but common mice such as the Microsoft Intellimouse and Microsoft Wheel Mouse will work just fine under Windows.

## Basic Mouse Control

On mice with a scroll wheel, you can push down on the wheel in order to use it as a middle button.

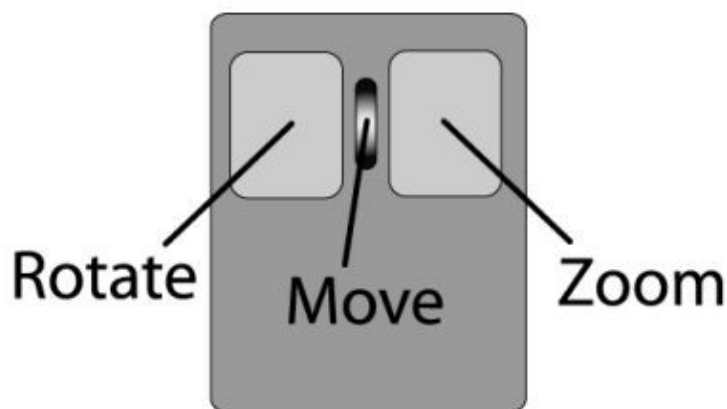

Here is a table of the basic mouse button/keyboard combinations for view manipulation:


Keyboard Modifier	Left Button	Middle Button	Right Button
(none)	Rotate Camera (Virtual Trackball)	Move Camera in XY (In Plane of Screen)	Move Camera in Z (Scale)
Shift Key			Move Clipping Planes
Control Keys			
Control and Shift Keys	Set Origin of Rotation		

An abbreviated version of this table, the Mouse Matrix, is always displayed in the Internal GUI, in order to help you remember which key and mouse button performs which action:

	L	M	R
None	Rota	Move	MovZ
Shft			Clip
Ctrl			
CtSh		Orig	

When using PyMOL on a laptop, it may be necessary to attach an external mouse or reassign the particular mouse controls you plan to use onto the reduced set of buttons that you have available internally (see reference on the "button" command).

## Virtual Trackball Rotation

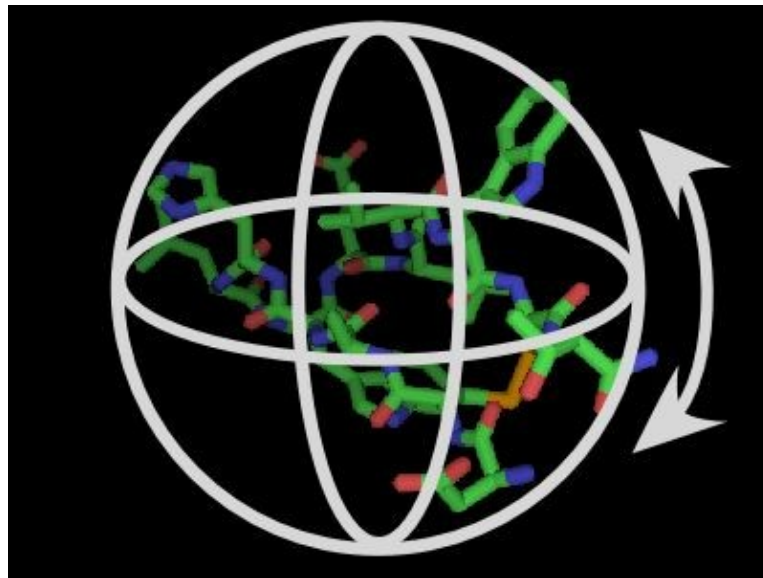


PyMOL's Virtual Trackball.

Virtual trackball rotation works as if there is an invisible ball in the center of the scene. When you click and drag on the screen, it is as if you put your finger on the sphere and rotated it in approximately the same manner. If you click outside the sphere, then you get rotation about the Z-axis only. Generally, the view will be easiest to control by either clicking in the center of the scene and moving outwards (mostly XY-rotation), or by clicking and dragging around the edge of the screen and moving in a circular fashion (Z-rotation).


## Moving Clipping Planes

PyMOL's clipping plane control is somewhat unusual and may take a few minutes to get used to. Instead of having separate controls for the front and back clipping planes, controls are combined into a single mode where **up–down mouse motion moves the front** (near) clipping plane and **left–right mouse motion controls the back** (far) clipping plane.

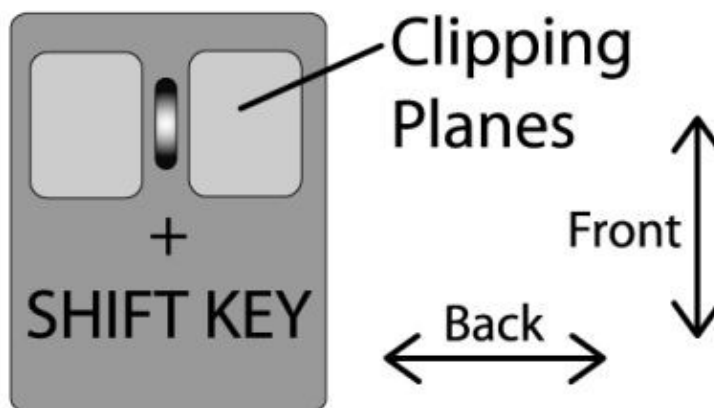

Control of clipping planes.

The advantage of the PyMOL clipping plane control is that tedious tandem manipulations of the clipping planes now becomes easy through the diagonal motions shown below.

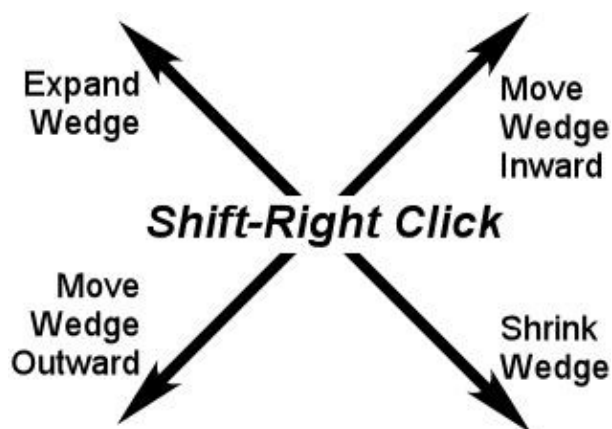


Changing the visible "wedge" by moving clipping planes in tandem.

You can also use the "clip" command to control the clipping planes.

## Changing the Origin of Rotation

When visualizing molecules, it is frequently necessary to change the origin of rotation so that you can inspect a particular region of the molecule. The fastest way to do this in PyMOL is to Control–Shift–Middle–Click on a visible atom in the scene.

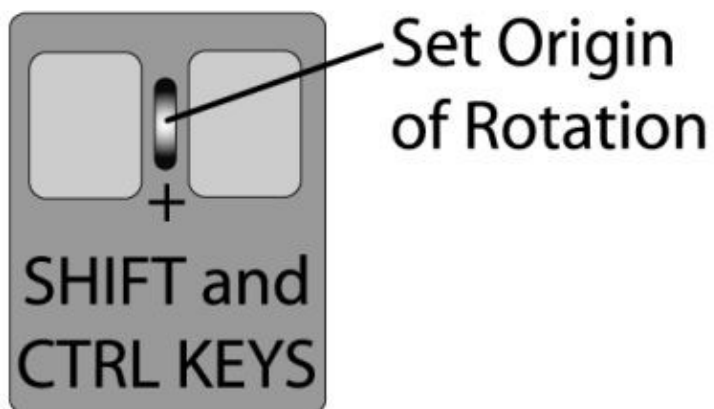

## Getting Comfortable

At this point, we recommend that you spend five or ten minutes getting comfortable with the controls described in this chapter. Specifically, you should be able to accomplish the following tasks:

1. Load a PDB file into PyMOL.
2. Rotate, translate, and zoom the camera.
3. Adjust the front and back clipping planes to clearly view a slice of the molecule.
4. Change the origin of rotation about any particular atom of interest.

# Getting Started with Commands

This section steps through a typical PyMOL session, introducing typed commands and describing how PyMOL responds to them. The details of command syntax are in the section titled "PyMOL Command Language."

The PyMOL language is case-sensitive, but upper case is not used in the current package. So just remember to type all your commands in lower case.

## Recording Your Work (Optional)

While you are learning PyMOL or doing complex projects, you may want to keep a record of all the commands you give in a plain text log-file that you can read and edit. To open a log-file, type the command **log\_open** followed by a file-name:

SYNTAX

```
log_open log-file-name
```

EXAMPLE

```
PyMOL> log_open log1.pml
```

All your commands, typed or clicked, will be recorded in the log-file. You should give your *log-file-name* the extension ".pml" so you can load the file as a script, to repeat your commands in a new session (see the subsection titled "Sessions and Scripts" below).

To stop recording your commands, type **log\_close**. If you don't type **log\_close** before you exit PyMOL, your log-file will still be saved to disk.

If you just want to save the current state of your PyMOL work without concern for the steps you took and the commands you gave, you can create a session-file (see "Sessions and Scripts").

## Loading Data

Next you need to input your data from a file, say atomic coordinates in PDB format:

SYNTAX

```
load data-file-name
```

EXAMPLE

```
PyMOL> load $PYMOL_PATH/test/dat/pept.pdb
```

Given this command, PyMOL will open and read the file "pept.pdb," create and name a corresponding object, display a representation of the object in the viewer, and add the object's name to the control panel.

By default, PyMOL names the object after the file it read. You can assign a different name to the object by typing the name in the command line:

SYNTAX

```
load data-file-name, object-name
```

#### EXAMPLES

```
PyMOL> load $PYMOL_PATH/test/dat/pept.pdb      # The object is named "pept".
                                                # PyMOL doesn't use
                                                # the file-name extension
                                                # ".pdb" in the object-name.

PyMOL> load $PYMOL_PATH/test/dat/pept.pdb, test # The object is named "test".
```

(Note that "#" is a comment character, so anything you type to the right of "#" in a command line is not interpreted by PyMOL.)

The command for loading a file follows typical PyMOL syntax. **load** is a *keyword*; it calls PyMOL to perform a certain function. *data-file-name* and *object-name* are arguments to **load**. These arguments tell PyMOL what file to load and what to name it, but in general, arguments to keywords just supply information that PyMOL needs to carry out commands.

## Manipulating Objects

After PyMOL creates an object, you can manipulate it in the view window and control panel with your mouse, and also by typing commands. For example, you can change from the default representation, called **lines**, to the more hefty **sticks**. First get rid of the **lines** and then show the **sticks**:

#### SYNTAX

```
hide representation
show representation
```

#### EXAMPLES

```
PyMOL> hide lines      # The object shown in lines disappears from view.

PyMOL> show sticks     # The object is represented as sticks in the viewer.
```

Other representations are **cartoons**, **ribbons**, **dots**, **spheres**, **surfaces**, and **meshes** (See the Section titled "Representations").

## Atom Selections

If you want to manipulate a subset of the atoms and bonds in a molecule, you can use *atom selections*. PyMOL is pretty sophisticated when it comes to selecting, grouping and naming groups of atoms. For example, you can select particular residues or atoms in a binding pocket, or a hydrophobic patch, or all the alanines in a helix, and so on. The Section titled "PyMOL Command Language" gives the details for selecting interesting groups of atoms.

You can use a selection just once, or you can name it to make it easier to use again later. For example, you can zoom in on a selection "on the fly:"

#### SYNTAX

```
zoom selection-expression  # Select the atoms just for zooming.
```

#### EXAMPLE

```
PyMOL> zoom resi 1-10      # The selector resi
                           # chooses amino acid residues
                           # given by the PDB sequence number
                           # identifier "1-10."
```

*Selection-expressions* range from single words to long complicated expressions. An *object-name* may be a *selection-expression*. The default *selection-expression* is **all**, which refers to all the atoms that are currently loaded. If a *selection-expression* is missing, PyMOL will apply the command to **all**. We'll keep our *selection-expressions* short in this section.

If you name the selection, you will be able to manipulate it any number of times. Object and selection names may include the upper or lower case characters (A/a to Z/z), numerals (0 to 9), and the underscore character (\_). Characters to avoid include:

```
! @ # $ % ^ & * ( ) ' " [ ] { } \ | ~ ` < > . ? /
```

First, name the selection:

#### SYNTAX

```
select selection-name, selection-expression
```

#### EXAMPLES

```
PyMOL> select akeeper, resi 1-10    # Select the residues and name them "akeeper."
```

Then use it:

#### SYNTAX

```
zoom selection-name
```

```
hide representation, selection-name
```

```
show representation, selection-name
```

#### EXAMPLES

```
PyMOL> zoom akeeper            # Zoom in on them in the viewer.
```

```
PyMOL> hide everything, akeeper # Hide their representation in the viewer.
```

```
PyMOL> show spheres, akeeper   # Show them in a different representation,
                                # spheres, this time.
```

When you create a *selection-name*, PyMOL puts it in the control panel so you can apply control panel functions to the selection using your mouse (See the section titled "PyMOL Command Language").

Named-selections such as "akeeper" are manipulated like PyMOL objects, but objects and named-selections are fundamentally different. PyMOL creates an *object-name* to locate data when you load a data file. Making selections is a way of pointing to a subset of that data. To distinguish *selection-names* from *object-names*, *selection-names* are shown inside parentheses in the control panel.



When you delete a *selection-name*, the data are still found under the *object-name*, but the data are no longer organized as the selection. In contrast, after you delete an object, you must reload the data to have access to it again.

#### SYNTAX

```
delete selection-name
```

```
delete object-name
```

#### EXAMPLES

```
PyMOL> delete akeeper           # "akeeper" is gone, but the object remains.
```

```
PyMOL> delete pept             # The atoms and bonds in "pept" are gone.
```

## Coloring Objects and Selections

You can apply various colors to selections and objects using typed commands. Predefined *color-names* are listed under the settings/colors pull-down menu. Many of them can be chosen from the control panel. See the section titled "Settings" to find out how to define more colors.

#### SYNTAX

```
color color-name                # All the representations of  
                                # loaded objects are colored.
```

```
color color-name, selection-expression  # The representation of  
                                         # the selection is colored.
```

#### EXAMPLES

```
PyMOL> color white              # Everything turns white.  
                                # This looks fine on the  
                                # default black background,  
                                # but causes disappearance  
                                # if you've changed the background to white.
```

```
PyMOL> color orange, pept       # Remember that "pept" is our object-name,  
                                # so the entire object is colored orange.
```

```
PyMOL> color green, resi 50+54+58  # The representation of  
                                # residues numbered 50, 54 and 58  
                                # is colored green.
```

```
PyMOL> color yellow, resi 60-90   # The representation of  
                                # residues numbered 60 through 90  
                                # is colored yellow.
```

```
PyMOL> color blue, akeeper        # Residues numbered 1-10,  
                                # which were collected in  
                                # the named selection "akeeper,"  
                                # are colored blue.
```

```
PyMOL> color red, ss h           # The representations of  
                                # helical residues  
                                # are colored red.
```

```
PyMOL> color yellow, ss s        # The representations of
```

```

# beta sheet residues
# are colored yellow.

PyMOL> color green, ss l+" "      # The representations of
# loop and unassigned residues
# are colored green.

```

In the last three examples, the *selector* **ss** chooses secondary structures specified by **h** for helix, **s** for beta sheet strand and **l+""** for loops and unspecified structures.

## Turning Objects and Selections On and Off

PyMOL can hold several objects in memory at the same time. The commands **disable** and **enable** allow you to eliminate representations of objects from the viewer while still controlling their properties with commands.

### SYNTAX

```
enable object-name
```

### EXAMPLE

```

PyMOL> load $PYMOL_PATH/test/dat/fc.pdb
PyMOL> load $PYMOL_PATH/test/dat/pept.pdb

PyMOL> disable pept              # All representations of "pept"
# are removed from view.

PyMOL> color yellow, name c+o+n+ca  # Backbone atoms in both "fc"
# and "pept" are colored yellow,
# but "pept" atoms
# are still not visible.

PyMOL> enable pept              # "pept" atoms are visible and
# its backbone atoms are yellow.

```

You can also use the **disable** command to get rid of the pink dots that identify the last-named selection in the viewer:

### SYNTAX

```
enable selection-name
```

### EXAMPLE

```

PyMOL> select bb, name c+o+n+ca    # Atoms included in the
# named-selection are indicated
# with pink dots in the viewer.

PyMOL> disable bb                 # The pink dots disappear,
# but the named selection "bb"
# is still visible.

PyMOL> color red, bb             # You can still manipulate "bb."

```

You can still operate on the representations of objects that are disabled, even with the commands **show** and **hide**. The results will be apparent when you subsequently **enable** the object.

## Changing Your Point of View

Dragging on a molecule with the mouse is often the easiest way to maneuver, but typed commands such as **zoom** and **orient** give you a different form of control, allowing computations to direct the view. **zoom**, as the name suggests, brings an object or selection close up in the center of the field of view. If the object or selection doesn't fit in the current view, the view opens out to include it. If it is just a small part of the current view, the view closes in to fill more of the screen with it.

### SYNTAX

```
zoom selection-expression      # The "camera" moves close
                               # to the selection so it fills the viewer,
                               # or moves further away to include
                               # all of the selection in the viewer
```

**orient** is a useful command when you want a fresh start in viewing the molecule. It aligns the object or selection so its largest dimension is shown horizontally, and its second largest dimension is shown vertically.

### SYNTAX

```
orient selection-expression    # The selection is aligned
                               # for maximum visibility in the viewer.
```

You can store orientations and recall them later in your PyMOL session using the command **view**. Storing a view only saves the viewpoint on the objects in the viewer. It does not save their representation. To store a view for later in the session, you need to "key" it, that is, to give a name or number as an argument to the command **view**. A second argument tells PyMOL whether you want it to store the view or recall it.

### SYNTAX

```
view key, action               # The possible actions are store and recall.
```

### EXAMPLES

```
PyMOL> view v1, store         # The current view is named "v1" and stored.

PyMOL> view v1, recall        # The view keyed "v1" is restored.

PyMOL> view v1                # recall is the default argument to view,
                               # so this also recalls "v1."
```

The keyword **view** only stores an orientation for the duration of the current PyMOL session. The next section gives the recipe for saving and restoring views in different PyMOL sessions.

## Saving Your Work

PyMOL saves your work in f kinds of processes: (1) Before you give a series of commands, you can initiate the process of logging your commands into plain text log-files that can later be used as scripts. (2) At any point in a PyMOL session you can save the memory state of the program by creating a session file that can later be loaded to restore that memory state. (3) You can write a graphics file to store the image you have created in the viewer for sharing or publication.

## Scripts and Log Files

A PyMOL script is just a text file, such as a log-file, containing typed PyMOL commands separated by carriage returns. When a script is loaded into PyMOL the commands it contains are executed. PyMOL expects scripts to have ".pml" file-name extensions (this is not strictly required, but it is good practice).

You can use log-files as scripts, and you can create scripts in a text editor such as **emacs**, **jot**, or **notepad**. It's often useful to keep a text editor open in a separate window while using PyMOL. Commands can then be cut and pasted between the two programs.

You can open a new log-file by typing **log\_open** *log-file-name*, or by clicking on "log" under the "File" menu and naming the log-file in the dialog box. You can also append commands to an existing log-file by choosing "append" or "resume" in the "File" menu. When you "resume" rather than "append," the existing log-file is first loaded as a script, and then subsequent commands are written to it.

Once you have opened a log-file in any of these ways, PyMOL will write and save all your commands, whether they are typed or given by clicking on the buttons in the GUI.

However, to store the orientation of a molecule into a log-file, you need to give the command **get\_view** (type it or use the GUI button). You may find it convenient to **get\_view** several times in a PyMOL session, and then edit the log-file to select the most useful views.

Scripts can be executed in several ways. Under Windows, scripts can be run in a new PyMOL session by double clicking on the script's icon. Alternatively, you can run a script using the "File" menu's "Run" option. PyMOL also understands "@" as the typed command that loads a script:

### SYNTAX

```
@script-file-name
```

### EXAMPLE

```
PyMOL> @my_script.pml
```

You can also include the *script-file-name* when launching PyMOL from a command line:

### SYNTAX

```
pymol script-file-name
```

### EXAMPLE (Windows)

```
C:\> pymol.exe my_script.pml
```

### EXAMPLE (Unix)

```
unix> ./pymol.com my_script.pml
```

## png Files

Once you are satisfied with the representation and orientation of your molecule, you may want to save the image in a graphics file. Before you do that, you can improve the quality of the graphics by switching from PyMOL's fast default graphics engine, OpenGL, to its ray tracer. The ray tracer is slower, but produces higher

quality renderings for display and publication. Ray tracing shows how light is reflected and how shadows fall in a three-dimensional world. Ray tracing may take some minutes for a large complex object. The keyword **ray** calls PyMOL's raytracer to redraw and display the image in the view window (See the section titled "Ray Tracing" for more details).

To save an image to a file, raytraced or not, use the "Save Image" option in the "File" menu or type the **png** command:

#### SYNTAX

```
png file-name
```

#### EXAMPLE

```
PyMOL> png $PYMOL_PATH/pep      # The file-name extension ".jpeg" is
                                # added. The image file "pep.jpeg" is stored
                                # in a path below the PyMOL installation.
```

The PNG file format is directly readable by PowerPoint. It can be converted into other formats using a package like ImageMagick.

## Session Files

If you want to be able to return to the current state of PyMOL, then you can create a session-file (Choose "Save Session" in the "File" menu and respond to the dialog box by naming the file with a ".pse" file-name extension). This utility works like the "Save" utility in a word processing program. A PyMOL session-file is a symbolic record of the state of PyMOL's memory, including the the objects that have been loaded or created, the named-selections that have been created, and the display in the viewer.

When you open the saved session-file, PyMOL's memory returns to the state that was saved. Because a session-file represents a PyMOL memory state, opening one means that you are eliminating everything that you currently have in PyMOL's memory, and replacing it with the memory state from the session-file.

A session-file differs from a log-file or a script in a number of ways. You have to open a log-file before you give the commands you want to save, but a session-file can be created at any point. A session-file is invoked by choosing "open" under the file menu, while a log-file is "run" as a script. Also, you can't write or edit session-files, as you can log-files and scripts.

It's a good idea to create session-files at strategic points in PyMOL sessions, for example, when you decide to explore one of several options. In this way, session-files can be used to replace an "undo" utility, which PyMOL lacks. You can store any number of PyMOL states in successive session-files, and revert to them, effectively "undo"-ing the work you did since creating the session-file.

## Command-Line Shortcuts

Since almost nobody likes to type, PyMOL's command-line interface includes several "shortcut" features designed to reduce typing. If you are a unix user, you will recognize the similarity with features found in tcsh or bash.

## Command Completion using TAB

If you type the first few characters of a command and then hit TAB, PyMOL will either complete the command or print out a list of which commands match the command.

EXAMPLE

```
PyMOL> sel
# hitting TAB will produce
PyMOL> select
```

If you hit the TAB key on a blank command line, PyMOL will output a list of its commands.

## Filename Completion using TAB

Some of the files you need to load into PyMOL may have long paths and filenames. PyMOL makes it easier to load such files by automatically completing unambiguous paths and filenames when you hit the TAB key. For instance,

EXAMPLE

```
PyMOL> load cry
# If "crystal.pdb" exists in the current directory, hitting TAB will generate
PyMOL> load crystal.pdb
```

If there is some ambiguity in the filename, PyMOL will complete the name up to the point of ambiguity and then print out the matching files in the directory.

## Automatic Inferences

There are a small number of "fixed string" arguments to PyMOL commands. For example, in

```
PyMOL> show sticks
```

"sticks" is a fixed string argument to **show**. Because there is only a small set of such arguments to **show**, PyMOL will infer your meaning even if you only provide it with a few letters. For example

```
PyMOL> show st
```

works just as well.

Keywords are also inferred in this manner, so

```
PyMOL> sh st
```

works too, as long as **show** is the PyMOL only command starting with "sh".

NOTE: PyMOL's command language continues to grow and develop, so it is important to use full-length commands and string arguments in scripts. Otherwise, you could not be sure that a later command or argument would not cause your abbreviation to become ambiguous. For example, "sh st" would no longer work if a **shutoff** command were added to the PyMOL language.

## Other Typed Commands and Help

This "Getting Started" section used the most frequent PyMOL commands in very brief examples. The section titled "Simple Examples" shows other commands that combine representations, selections and property changes. More complicated examples appear in the section titled "Cookbook and Complex Examples," and a comprehensive listing of typed commands appears in the section titled "Command and API Reference."

To see a list of the keyword commands that are available in PyMOL on your computer screen, type **help** and "enter" (Typing TAB and "enter" will work too). Add the keyword if you want help on a particular command:

SYNTAX

```
help keyword
```

EXAMPLE

```
PyMOL> help load
```

PyMOL responds by displaying the manual page that describes the command in the PyMOL viewer. Command line completion works inside of help, so if you don't remember the full keyword, type **help**, the first character or so of the keyword, and hit TAB. Python will display a list of possible help topics.

Click inside the viewer and hit escape to toggle back and forth between the display of the manual page or the list of commands and the molecules you have loaded in PyMOL.

All the keywords that PyMOL understands are listed alphabetically and described in the "Reference" section. PyMOL commands run on top of the Python programming language and may contain Python statements. After you type in a command and hit return, PyMOL will check whether the first word is one of its keywords (or if it can be extended into a keyword). If not, PyMOL will pass the command on to the Python interpreter. PyMOL will return a Python error message if neither a PyMOL nor a Python keyword is recognized.

# Command Syntax and Atom Selections

## Syntax

A typed PyMOL command always starts with a keyword that calls PyMOL to execute an action. It ends with a carriage return ("enter" on your keyboard).

The simplest commands consist of a keyword alone. For example, typing **quit** will end your PyMOL session. The **quit** command never takes an argument.

Many commands have default arguments, so you can type the keyword alone and PyMOL will supply the rest. For example, the default argument to **zoom** is the *selection-expression* **all**:

EXAMPLE

```
PyMOL> zoom                # All visible representations
                             # are included in the view.
```

For some keywords, certain arguments are required and others are supplied by default. For example, the keyword **color** requires one argument, the *color-name*. As for **zoom**, the default *selection-expression* is **all**:

SYNTAX

```
color color-name

color color-name, selection-expression
```

EXAMPLES

```
PyMOL> color red            # All the representations
                             # are colored red.

PyMOL> color red, name ca    # Only the representations of
                             # atoms named c-alpha are colored red.
```

When you type a command that has more than one argument, **color** *color-name*, *selection-expression* in this case, a comma must separate the arguments.

*Selection-expressions* are an essential type of keyword argument. They can be simple or complex, with several different kinds of syntax.

## Selection-expressions

*Selection-expressions* stand for lists of atoms in arguments that are subject to PyMOL commands. You can name the selections to facilitate their re-use, or you can specify them anonymously (without names). Object and selection names may include the upper or lower case characters A/a to Z/z, numerals 0 to 9, and the underscore character (\_). Characters to avoid include:

```
! @ # $ % ^ & * ( ) ' " [ ] { } \ | ~ ` < > . ? /
```

*Selection-expressions* describe the class of atoms you are referencing. Most of them require identifiers to complete the specification. For example, the selector **resi** references biopolymer residues by sequence



number, and the identifier gives the number. The selector **name** references atoms according to the names described in the PDB, and the identifier gives the name (**ca** for alpha carbons, **cb** for beta carbons, etc). A handful of *selection-expressions* don't require identifiers, but most do.

PyMOL uses several logical operators to increase the generality or specificity of *selection-expressions*. Logical combinations of selectors can get complex, so PyMOL accepts short forms and macros that express them with a minimum of keystrokes. This section describes named-selections, and then gives the syntax for making selections in a progression from simple one-word selectors to complex combinations of selectors, using macros and short forms.

## Named Atom Selections

Atom selections can be named for repeated use by using the **select** command:

### SYNTAX

```
select  selection-name, selection-expression
                                     # The  selection-name and
                                     # the  selection-expression
                                     # are both arguments to select
                                     # so they are separated by a comma.
```

### EXAMPLE

```
PyMOL> select bb, name c+o+n+ca    # Create an atom selection named "bb"
                                     # including all atoms named
                                     # "C","O","N", or "CA";
PyMOL> color red, bb                # color the selection red,
PyMOL> hide lines, bb               # hide the line representation,
PyMOL> show sticks, bb              # show it using the stick representation,
PyMOL> zoom bb                      # and zoom in on it.
```

In this case, the *selection-expression* is the property selector **name**, which takes the list of identifiers **ca+c+n+o** to complete the specification. Property selectors and their identifiers are discussed below.

Named atom selections appear in the PyMOL names list in the control panel. They are distinguished from objects by a surrounding set of parentheses. The control panel options available under the diamond menu differ between objects and atom-selections, because objects and named selections play slightly different roles in PyMOL. Named selections are pointers to subsets of data that are found under an object name. After an object is deleted, the data are no longer available, unless you reload the object. Any named selections that refer to atoms in that object will no longer work. But when named selections are deleted, the data are still available under the object name. Disabling objects eliminates them from the viewer, but disabling named-selections just turns off the pink dots that highlight them in the viewer.

Atom-selections, named or not, can span multiple objects:

### EXAMPLE

```
PyMOL> load $PYMOL_PATH/test/dat/fc.pdb
PyMOL> load $PYMOL_PATH/test/dat/pept.pdb

PyMOL> select alpha_c, name ca      # The named selection "alpha_c"
                                     # is created -- it includes atoms
                                     # in both "fc" and "pept" objects.
PyMOL> color red, name ca           # "CA" atoms in both objects
```

# are colored red.

Named selections will continue working after you have made changes to a molecular structure:

#### EXAMPLE

```
PyMOL> load $PYMOL_PATH/test/dat/pept.pdb
PyMOL> select bb, name c+o+n+ca      # The named selection "bb"
                                      # is created.

PyMOL> count_atoms bb                # PyMOL counts 52 atoms in "bb."

PyMOL> remove resi 5                 # All atoms in residue 5 are removed
                                      # from the object "pept."

PyMOL> count_atoms bb                # Now PyMOL counts
                                      # the remaining 48 atoms in "bb."
```

Named selections are static. Only atoms that exist at the time the selection is defined are included in the selection, even if atoms which are loaded subsequently fall within the selection criterion:

#### EXAMPLE

```
PyMOL> load $PYMOL_PATH/test/dat/pept.pdb

PyMOL> select static_demo, pept      # The named selection "static_demo"
                                      # is created to reference all atoms.

PyMOL> count_atoms static_demo        # PyMOL counts 107 atoms
                                      # in "static_demo."

PyMOL> h_add                          # PyMol adds hydrogens in
                                      # the appropriate places

PyMOL> count_atoms static_demo        # PyMOL still counts 107 atoms
                                      # in "static_demo,"
PyMOL> count_atoms                    # even though it counts 200 atoms
                                      # in "pept."
```

Named selections can also be used in subsequent atom selections:

#### EXAMPLE

```
PyMOL> select bb, name c+o+n+ca      # An atom selection named "bb"
                                      # is made, consisting of all
                                      # atoms named "C","O","N", or "CA."

PyMOL> select c_beta_bb, bb or name cb
                                      # An atom selection named "c_beta_bb"
                                      # is made, consisting of
                                      # all atoms named "C", "O", "N", "CA" or "CB."
```

Note that the word "or" is used to select all atoms in the two groups, "bb" and "cb." The word "and" would have selected no atoms because it is interpreted in its boolean logical sense, not its natural language sense. See the subsection on "Selection Algebra" below.

## Single-word Selectors

The very simplest *selection-expressions* are single-word selectors. These selectors do not take identifiers; they are complete by themselves.

Single Word Selector	Short Form Selector	Description
all	*	All atoms currently loaded into PyMOL
none	none	No atoms (empty selection)
hydro	h.	All hydrogen atoms currently loaded into PyMOL
hetatm	het	All atoms loaded from Protein Data Bank HETATM records
visible	v.	All atoms in enabled objects with at least one visible representation
present	pr.	All atoms with defined coordinates in the current state (used in creating movies)

The selector **none** won't come up much when you are typing commands directly into PyMOL, but it is useful in programming scripts.

As the table shows, many single-word selectors have short forms to save on typing. Some short forms must be followed by a period and a space, in order to delimit the word. Short forms and long forms have the same effect, so choose the form that suits you:

### EXAMPLES

```
PyMOL> color blue, all           # It all turns blue.
PyMOL> color blue, *

PyMOL> hide hydro                # Representations of all
PyMOL> hide h.                  # hydrogen atoms are hidden.

PyMOL> show spheres, hetatom    # All the atoms defined as HETATOMS
PyMOL> show spheres, het        # in the PDB input file
                                   # are represented as spheres.
```

## Property Selectors

PyMOL reads data files written in PDB, MOL/SDF, Macromodel, ChemPy Model, and Tinker XYZ formats. Some of the data fields in these formats allow PyMOL to assign properties to atoms. You can group and select atoms according to these properties using property selectors and identifiers: the selectors correspond to the fields in the data files, and the identifiers correspond to the target words to match, or the target numbers to compare.

The items in a list of identifiers are separated by plus signs (+) only. Do not add spaces within a list of identifiers. The selector **resi** takes (+)-separated lists of identifiers, as in

### EXAMPLE

```
PyMOL> select nterm, resi 1+2+3
```

or, alternatively, it may take a range given with a dash:

### EXAMPLE

```
PyMOL> select nterm, resi 1-3
```

However, you will get an error message if you try to combine a list and a range in an identifier to a **resi** as in "select mistake, resi 1-3+6."

The identifier for a blank field in an input file is an empty pair of quotes:

EXAMPLE

```
PyMOL> select unstruct, ss ""      # A named selection is created
                                     # to contain all atoms that are not assigned
                                     # a secondary structure.
```

Most *property selectors* select matches to their identifiers:

Matching Property Selector	Short Form Selector	Identifier and Example
symbol	e.	<i>chemical-symbol-list</i> list of 1- or 2-letter chemical symbols from the periodic table  PyMOL> select polar, symbol o+n
name	n.	<i>atom-name-list</i> list of up to 4-letter codes for atoms in proteins or nucleic acids  PyMOL> select carbons, name ca+cb+cg+cd
resn	r.	<i>residue-name-list</i> list of 3-letter codes for amino acids  PyMOL> select aas, resn asp+glu+asn+gln  or list of up to 2-letter codes for nucleic acids  PyMOL> select bases, resn a+g
resi	i.	<i>residue-identifier-list</i> list of up to 4-digit residue numbers  PyMOL> select mults10, resi 1+10+100+1000  <i>residue-identifier-range</i>  PyMOL> select nterm, resi 1-10
alt	alt	<i>alternate-conformation-identifier-list</i> list of single letters  PyMOL> select altconf, alt a+""
chain	c.	<i>chain-identifier-list</i> list of single letters or sometimes numbers

		<b>PyMOL&gt;</b> select firstch, chain a
segi	s.	<i>segment-identifier-list</i> list of up to 4 letter identifiers  <b>PyMOL&gt;</b> select ligand, segi lig
flag	f.	<i>flag-number</i> a single integer from 0 to 31  <b>PyMOL&gt;</b> select f1, flag 0
numeric_type	nt.	<i>type-number</i> a single integer  <b>PyMOL&gt;</b> select type1, nt. 5
text_type	tt.	<i>type-string</i> a list of up to 4 letter codes  <b>PyMOL&gt;</b> select subset, text_type HA+HC
id	id	<i>external-index-number</i> a single integer  <b>PyMOL&gt;</b> select idno, id 23
index	idx.	<i>internal-index-number</i> a single integer  <b>PyMOL&gt;</b> select intid, index 11
ss	ss	<i>secondary-structure-type</i> list of single letters  <b>PyMOL&gt;</b> select allstrs, ss h+s+l+" "


Other *property selectors* select by comparison to numeric identifiers:

<b>Numeric Selector</b>	<b>Short Form</b>	<b>Argument &amp; Example</b>
b	b	<i>comparison-operator b-factor-value</i> a real number  <b>PyMOL&gt;</b> select fuzzy, b > 10
q	q	<i>comparison-operator occupancy-value</i> a real number  <b>PyMOL&gt;</b> select lowcharges, q <0.50
formal_charge	fc.	<i>comparison-operator formal charge-value</i> an integer  <b>PyMOL&gt;</b> select doubles, fc. = -1
partial_charge	pc.	<i>comparison-operator partial charge-value</i> a real number  <b>PyMOL&gt;</b> select hicharges, pc. > 1

Details of the atom and residue name formats can be found in the official guide to PDB file formats, [http://www.rcsb.org/pdb/docs/format/pdbguide2.2/guide2.2\\_frame.html](http://www.rcsb.org/pdb/docs/format/pdbguide2.2/guide2.2_frame.html).

## Selection Algebra

Selections can be made more precise or inclusive by combining them with logical operators, including the boolean **and**, **or** and **not**. The boolean **and** selects only those items that have both (or all) of the named properties, and the boolean **or** selects items that have either (or any) of them. Venn diagrams show that **and** selects the areas of overlap, while **or** selects both areas.

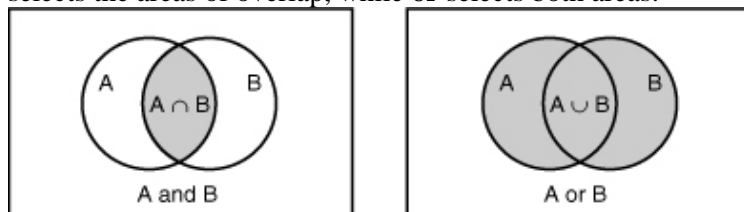

### Operators:

Selection operators and modifiers are listed below. The dummy variables *s1* and *s2* stand for selection-expressions such as "chain a" or "hydro."

Operator	Short form	Effect
not <i>s1</i>	! <i>s1</i>	Selects atoms that are not included in <i>s1</i>  <code>PyMOL&gt; select sidechains, ! bb</code>
<i>s1</i> and <i>s2</i>	<i>s1</i> & <i>s2</i>	Selects atoms included in both <i>s1</i> and <i>s2</i>  <code>PyMOL&gt; select far_bb, bb &amp; farfm_ten</code>
<i>s1</i> or <i>s2</i>	<i>s1</i>   <i>s2</i>	Selects atoms included in either <i>s1</i> or <i>s2</i>  <code>PyMOL&gt; select all_prot, bb   sidechain</code>
<i>s1</i> in <i>s2</i>	<i>s1</i> in <i>s2</i>	Selects atoms in <i>s1</i> whose identifiers name, resi, resn, chain and segi all match atoms in <i>s2</i>  <code>PyMOL&gt; select same_atms, pept in prot</code>
<i>s1</i> like <i>s2</i>	<i>s1</i> l. <i>s2</i>	Selects atoms in <i>s1</i> whose identifiers name and resi match atoms in <i>s2</i>  <code>PyMOL&gt; select similar_atms, pept like prot</code>
<i>s1</i> gap <i>X</i>	<i>s1</i> gap <i>X</i>	Selects all atoms whose van der Waals radii are separated from the van der Waals radii of <i>s1</i> by a minimum of <i>X</i> Angstroms.  <code>PyMOL&gt; select farfm_ten, resi 10 gap 5</code>
<i>s1</i> around <i>X</i>	<i>s1</i> a. <i>X</i>	Selects atoms with centers within <i>X</i> Angstroms of the center of any atom in <i>s1</i>  <code>PyMOL&gt; select near_ten, resi 10 around 5</code>
<i>s1</i> expand <i>X</i>	<i>s1</i> e. <i>X</i>	Expands <i>s1</i> by all atoms within <i>X</i> Angstroms of the center of any atom in <i>s1</i>  <code>PyMOL&gt; select near_ten_x, near10 expand 3</code>

<i>s1</i> within X of <i>s2</i>	<i>s1</i> w. X of <i>s2</i>	Selects atoms in <i>s1</i> that are within X Angstroms of the <i>s2</i>  <b>PyMOL&gt;</b> select bbnearten, bb w. 4 of resi 10
byres <i>s1</i>	br. <i>s1</i>	Expands selection to complete residues  <b>PyMOL&gt;</b> select complete_res, br. bbnear10
byobject <i>s1</i>	bo. <i>s1</i>	Expands selection to complete objects  <b>PyMOL&gt;</b> select near_obj, bo. near_res
neighbor <i>s1</i>	nbr. <i>s1</i>	Selects atoms directly bonded to <i>s1</i>  <b>PyMOL&gt;</b> select vicinos, neighbor resi 10

Logical selections can be combined. For example, you might select atoms that are part of chain a, but not residue number 125:

#### EXAMPLE

```
PyMOL> select chain a and (not resi 125)      # selects atoms that are part of
                                              # chain a, but not
                                              # residue number 125.

PyMOL> select (name cb or name cg1 or name cg2) and chain A  # These two
                                                              # selections are
PyMOL> select name cb+cg1+cg2 and chain A      # equivalent.
                                                              # select c-beta's,
                                                              # c-gamma-1's and
                                                              # c-gamma-2's
                                                              # that are
                                                              # in chain A.
```

Like the results of groups of arithmetic operations, the results of groups of logical operations depend on which operation is performed first. They have an order of precedence. To ensure that the operations are performed in the order you have in mind, use parentheses:

```
byres ((chain a or (chain b and (not resi 125))) around 5)
```

PyMOL will expand its logical selection out from the innermost parentheses.

## Atom Selection Macros

Macros make it possible to represent a long atom selection phrase such as

```
PyMOL> select pept and segi lig and chain b and resi 142 and name ca
```

in a more compact form:

```
PyMOL> select /pept/lig/b/142/ca
```

An atom selection macro uses slashes to define fields corresponding to identifiers. The macro is used to select atoms using the boolean "and," that is, the selected atoms must have all the matching identifiers:

```
/object-name/segi-identifier/chain-identifier/resi-identifier/name-identifier
```

These identifiers form a hierarchy from the object-name at the top, down to the name-identifier at the bottom. PyMOL has to be able to recognize the macro as one word, so no spaces are allowed within it.

Macros come in two flavors: those that begin with a slash and those that don't. The presence or absence of a slash at the beginning of the macro determines how it is interpreted. If the macro begins with a slash, PyMOL expects to find the fields *starting* from the *top* of the hierarchy: the first field to the right of the slash is interpreted as an object-name; the second field as an identifier to segi; the third as an identifier to chain, and so on. It may take any of the following forms:

```
/object-name/segi-identifier/chain-identifier/resi-identifier/name-identifier
/object-name/segi-identifier/chain-identifier/resi-identifier
/object-name/segi-identifier/chain-identifier
/object-name/segi-identifier
/object-name
```

#### EXAMPLES

```
PyMOL> zoom /pept
PyMOL> show spheres, /pept/lig/
PyMOL> show cartoon, /pept/lig/a
PyMOL> color pink, /pept/lig/a/10
PyMOL> color yellow, /pept/lig/a/10/ca
```

If the macro does not begin with a slash, it is interpreted differently. In this case, PyMOL expects to find the fields *ending* with the *bottom* of the hierarchy. Macros that don't start with a slash may take the following forms:

```
resi-identifier/name-identifier
chain-identifier/resi-identifier/name-identifier
segi-identifier/chain-identifier/resi-identifier/name-identifier
object-name/segi-identifier/chain-identifier/resi-identifier/name-identifier
```

#### EXAMPLES

```
PyMOL> zoom 10/cb
PyMOL> show spheres, a/10-12/ca
PyMOL> show cartoon, lig/b/6+8/c+o
PyMOL> color pink, pept/enz/c/3/n
```

You can also omit fields between slashes. Omitted fields will be interpreted as wildcards, as in the following forms:

```
resi-identifier/
resi-identifier/name-identifier
chain-identifier//
object-name//chain-identifier
```

#### EXAMPLES

```
PyMOL> zoom 142/                # Residue 142 fills the viewer.

PyMOL> show spheres, 156/ca      # The alpha carbon of residue 156
                                # is shown as a sphere

PyMOL> show cartoon, a//         # Chain "A" is shown as a cartoon.

PyMOL> color pink, pept//b       # Chain "B" in object "pept"
                                # is colored pink.
```



Selection macros must contain at least one forward slash in order to distinguish them from other words in the selection language. Being words, they must not contain any spaces. When using macros, it is also important to understand that they are converted into long form before being submitted to the selection engine. This can help in the interpretation of error messages.

## Calling Python from within PyMOL

Single-line Python statements can be issued directly within PyMOL. For example:

```
PyMOL> print 1 + 2
3
```

Full access is available to standard Python library functions, and you can assign results to symbols.

```
PyMOL>import time
PyMOL>now = time.time()
PyMOL>print now
1052982734.94
```

Multi-line blocks of Python can be included within PyMOL command scripts provided that a backslash ('\') is used for continuation on all but the final line.

```
PyMOL> for a in range(1,6): \
PyMOL>     b = 6 - a \
PyMOL>     print a, b
1 5
2 4
3 3
4 2
5 1
```

# Cartoon Representations

## Background

### Accessibility

Cartoon ribbons in PyMOL rival those of the popular Molscript/Raster3D packages, but PyMOL makes creating high quality images much easier. While PyMOL can read Molscript output directly (see the chapter on Molscript), this is no longer necessary or as convenient as utilizing PyMOL's built-in cartoon ribbon capability:



PyMOL built-in ribbons




"molauto -nice ... | molscript -r > ..."

Molscript's cartoons are slightly more ideal, but PyMOL comes pretty darn close!

Note that all of the images in this section were colored using the rainbow feature (Color pop-up menu) and ray-traced with antialiasing enabled.

### Pretty *and* Correct


One of the advantages of PyMOL's cartoon ribbon facility is that it is easy to switch between "smoothed" versions of protein secondary structure, and "correct" renditions which portray actual main chain coordinates. Although cartoons are often used solely to represent protein structures in a schematic sense, sometimes it is desirable to combine a schematic overall picture with atomic resolution in one particular location. However, unless the cartoon track properly with alpha-carbon positions, the resulting figures will look a little silly:

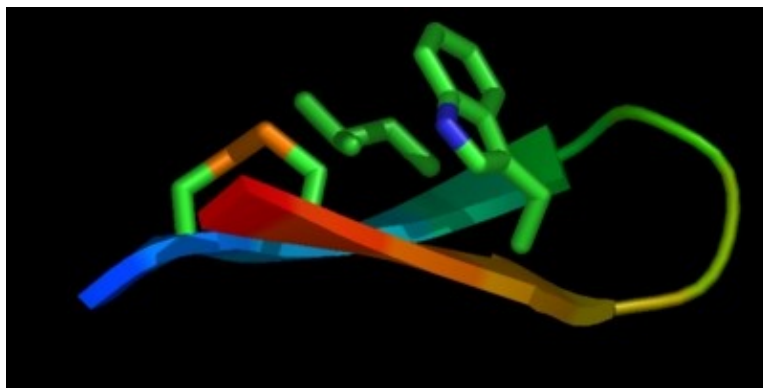

In the above image, the side chains are floating off into space. Disabling "flat sheets" from the Cartoons Menu or issuing the command


```
set cartoon_flat_sheets, 0
```

will make the beta strands follow the true path of the backbone through space and give a more accurate rendition of the structure.

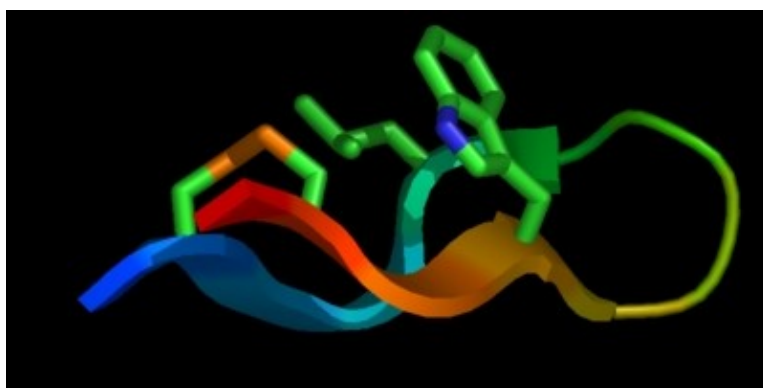


The appearance of a cartoon over the entire molecule will be substantially different when all smoothing features are turned off. For instance, with smoothing enabled:

```
set cartoon_flat_sheets, 1
set cartoon_smooth_loops, 1
```

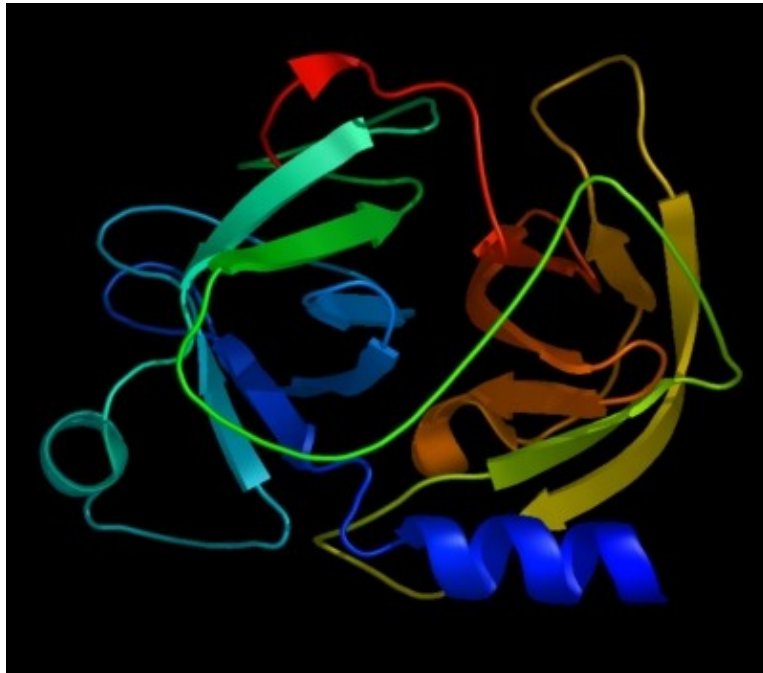


the image differs substantially from:

```
set cartoon_flat_sheets, 0  
set cartoon_smooth_loops, 0
```

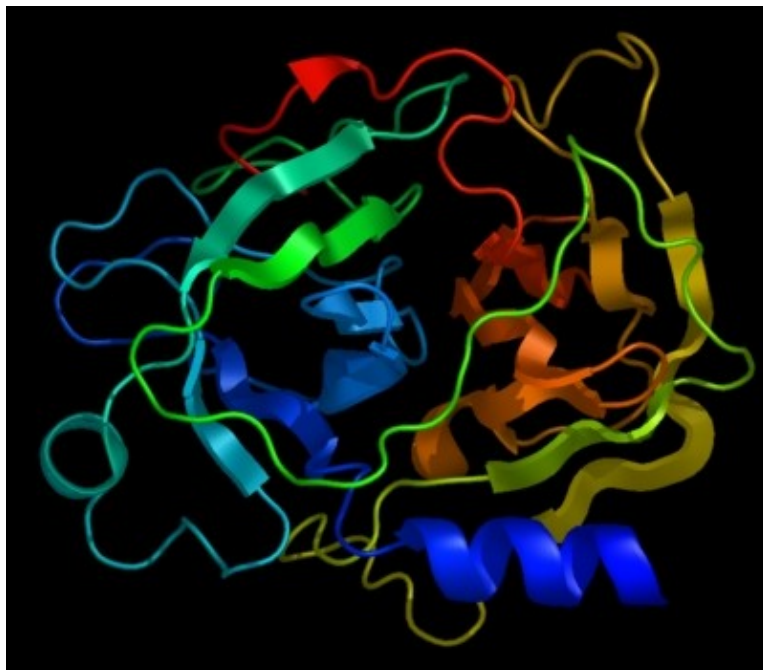

which more accurately reflects the true path of the peptide backbone:

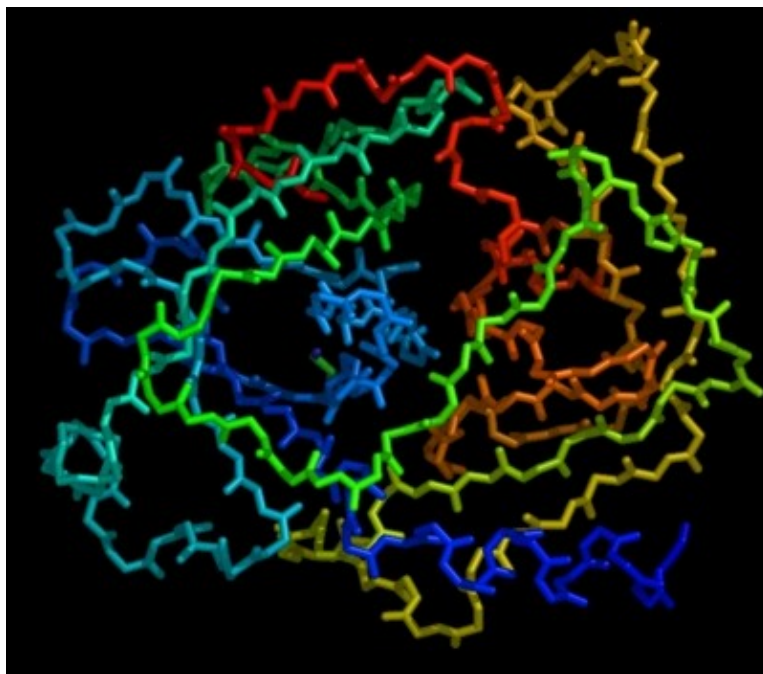

To facilitate beautiful imagery, smoothing is enabled by default (just like Molscript) [NOTE: THIS MAY CHANGE BEFORE VERSION 1.0] . Just be sure to turn it off when you want to study structures at atomic resolution (*remember, real life is a bit more complicated than what you see in cartoons!*).

## Customization

### Cartoon Types


Best results will be obtained when secondary structure information has been defined for each residue in the molecule. Under these conditions, PyMOL will do extra processing to insure that good normals have been calculated for helical regions, and perform smoothing of loops, where desired.

Also under such conditions, in automatic mode, cartoon representations will be assigned according to the secondary structure type. However, you can instruct PyMOL to ignore such information, and manually control when and where various cartoon representations are employed.


```
show cartoon
cartoon automatic      # default
```




cartoon loop

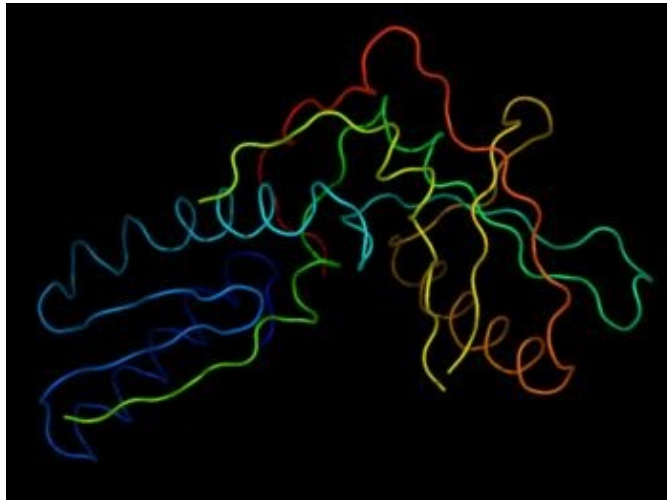


cartoon rect

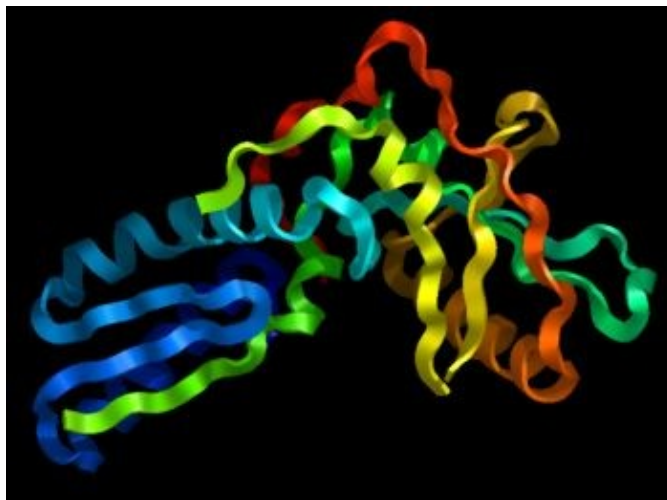

cartoon oval

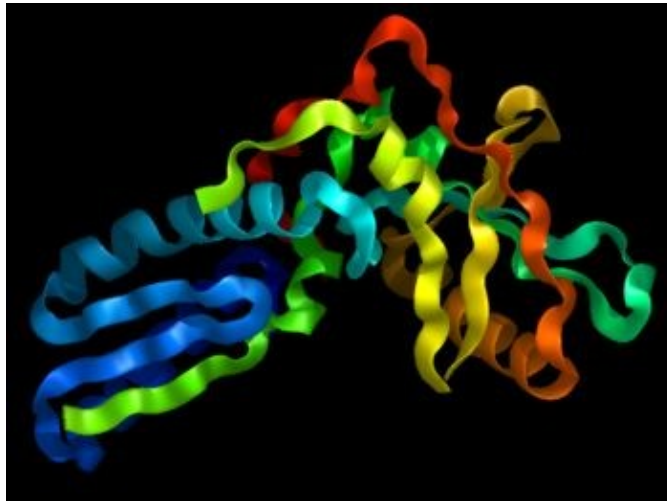


cartoon tube

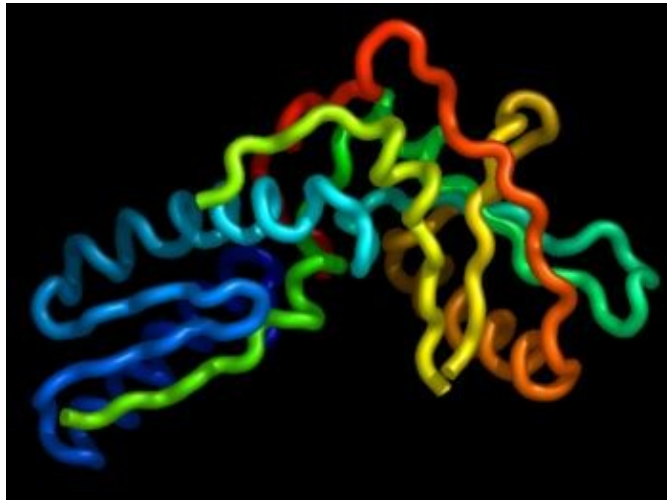

```
cartoon tube, 1:49/  
cartoon arrow, 50:99/  
cartoon loop, 100:149/  
cartoon oval, 150:199/  
cartoon rect, 200:250/
```




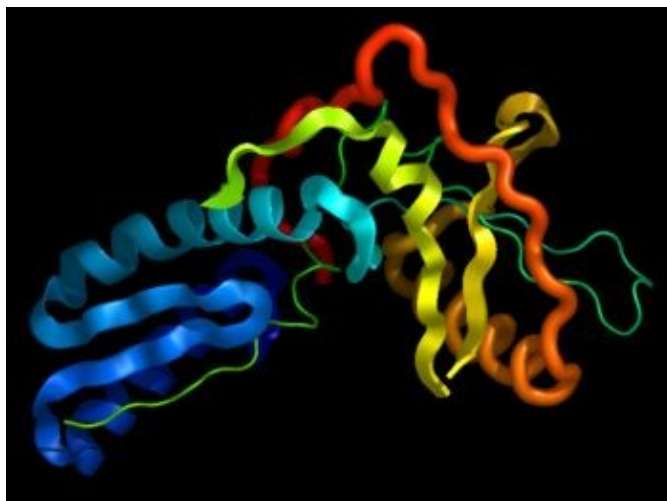

All cartoon ribbons have associated parameters accessible from the "set" command which allow you to change their appearance. See the chapter on Settings for more information.

## Fancy Helices


```
set cartoon_fancy_helices, 0
```

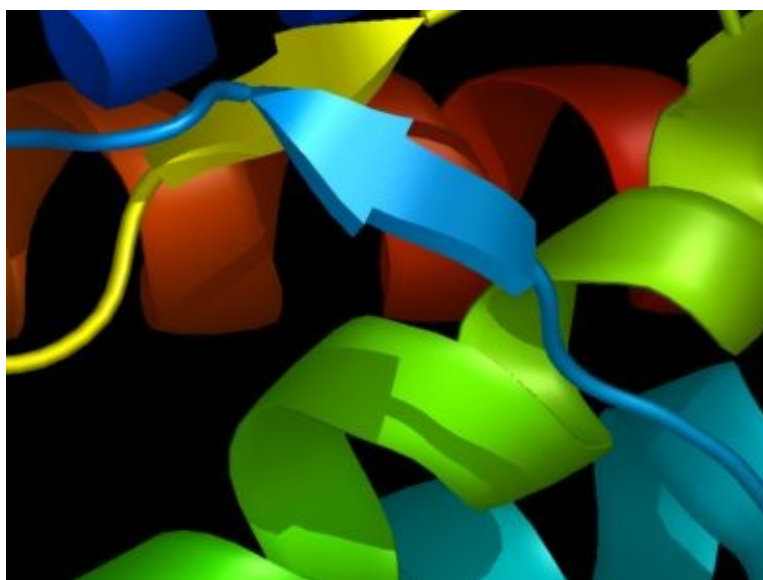

Molscript addicts who simply must have those ribbon helices with tubular edges will not be disappointed with "fancy helices":

```
set cartoon_fancy_helices, 1
```



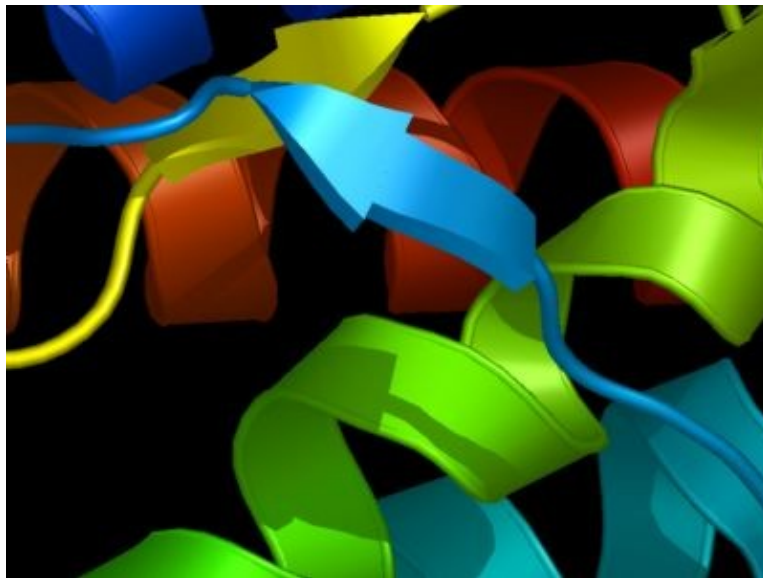

## Secondary Structure Assignment

It is recommended that you read in PDB files which already have correct secondary structure assignments from a program like DSSP. However, PyMOL does have a reasonably fast secondary structure alignment algorithm called "dss". Please be aware that due to the subjective nature of secondary structure assignment in borderline cases, dss results will differ somewhat from DSSP.

### SYNTAX

```
dss selection
```

### EXAMPLE

```
dss ldfr
```

If you are visualizing an animation, you may wish derive secondary structure assignment from a specific state of the animation. This can be done with:

### SYNTAX

```
dss selection, state
```

### EXAMPLE


```
dss mov, 1
```

To change assignments manually, the best way is to use the alter command as follows:

```
show cartoon
alter 11-40/, ss='H'           # assign residues 11-40 as helix
alter 40-52/, ss='L'           # assign residues 40-52 as loop
alter 52-65/, ss='S'           # assign residues 52-65 as sheet
alter 65-72/, ss='H'           # assign residues 65-72 as helix
rebuild                       # regenerate the cartoon
```

# Ray-Tracing

Ray-tracing produces the highest quality molecular graphics images. **PyMOL is the first full-featured molecular graphics program to include a high-speed ray-tracer** which works with its native internal geometries (except text).

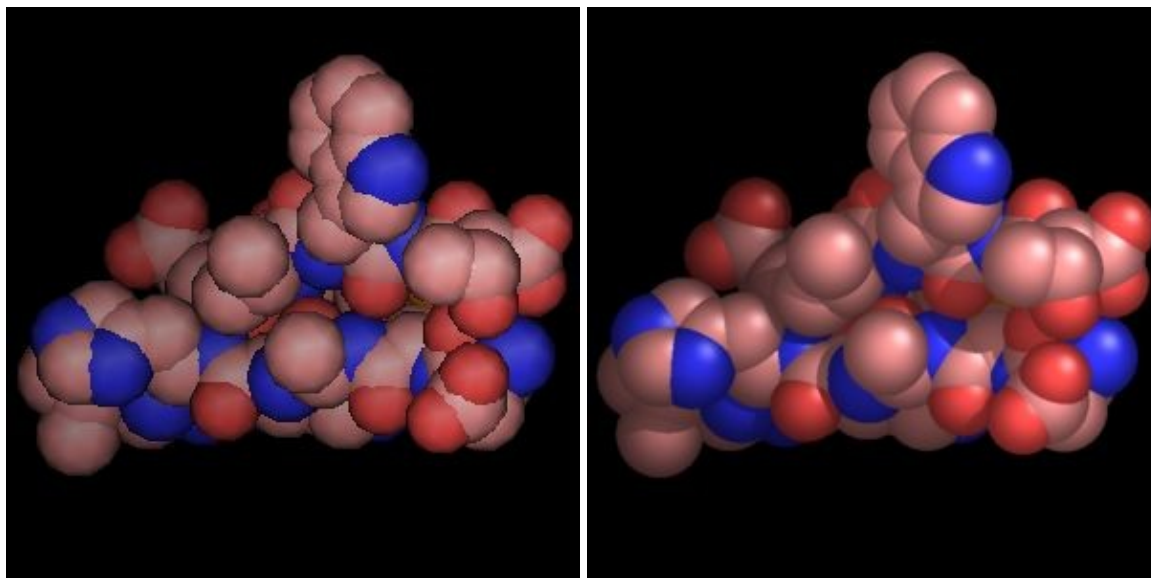

OpenGL Rendering (real-time manipulation)      Ray-traced Rendering (seconds or minutes per frame)  
You can ray-trace any Scene in PyMOL by clicking the "Ray Trace" button in the external GUI or using the "ray" command. The built-in raytracer also makes it possible to easily assemble very high-quality movies in a snap.

## Important Settings

These can be changed using the "set" command. Unless otherwise specified, the settings apply only to the ray-tracing engine and not the OpenGL renderer. Some reconciliation between the two renderers is much needed, so be warned that these settings may change in the future.

Normally, the only settings you will need to change are **orthoscopic**, **antialias**, and **gamma**. If you are down in an enzyme active site which is heavily shadowed, you may want to increase **direct** to 0.5–0.7 in order to improve brightness and contrast.

- **orthoscopic** (0 or 1) controls whether the OpenGL renderer uses the same orthoscopic transformation as the renderer. You'll want to set this to 1 when preparing figures so that OpenGL and raytracing match pixel-for-pixel.
- **ambient** (0.0–1.0) controls the ambient light intensity for both OpenGL and the ray-tracer.
- **ambient\_scale** (float) controls the relative ambient intensity between OpenGL and the ray-tracer.
- **antialias** (0 or 1) generate a "smooth" image (best quality, but takes 4X as long).
- **direct** (0.0–1.0) the planer light intensity originating from the camera.

- **gamma** (0.1–2.0) gamma transformation applied after rendering is complete.
- **light** (vector) the position of the light.
- **reflect** (0.0–1.0) the planer light intensity originating from the light source.
- **spec\_reflect** (0.0–1.0) intensity of the specular reflection from the light.
- **spec\_power** (1–100) how localized is the specular reflection (higher=smaller).

## Saving Images

### png

All images (ray-traced or not) can be saved in PNG format using the "png" command. This format is directly readable by PowerPoint, and can be easily converted into other formats using a package like ImageMagick. You can also save images using the "Save Image" option in the "File" menu. Images are always saved at the same resolution as the viewer window.

```
ray
png my_image.png
```

# Stereo

## Introduction

PyMOL can supports several different stereo graphics options.

## Supported Stereo Modes

**Crosseye Stereo**

**Walleye Stereo**

**Hardware Stereo**

## Generating Stereo Figures

Stereo figures are often used in molecular graphics illustrations in order to provide the reader with a three-dimensional view using a two-dimensional sheet of paper. To achieve a satisfying 3D effect, care must be taken to insure that the stereo pairs are generated with the proper transformation and printed with the correct separation between images.

A key challenge in preparation of stereo figures is the fact that journals usually shrink images before they are printed. In preparing your figure, it is essential that your drafts correspond to the final printed size. You may need to obtain this information from the journal ahead of time.

# Movies

## Concepts

### States and Frames

PyMOL has a powerful and unique molecular movie-making capability. In order to use it, you first need to understand a few key concepts:

- **States:** States most directly correspond to particular arrangements of atoms at a point in time. For example, they could consist of steps in a molecular dynamics simulation or individual points of a coordinate interpolation. If you are making a movie of a static coordinate set (such as a single crystal structure) then you have only one state. All objects in PyMOL can potentially consist of multiple states.
- **Frames:** Frames are like the individual images you'd find on a movie reel, except that in PyMOL, frames are composed of states instead of images, and frames can have additional actions associated with them (such as rotation of the camera).

The user can fully interact with models while movies are playing.

**NOTE:** State and frame indexes begin with 1, and not 0 as most C and Python programmers would expect. If you load states into an object with a state index of 0, you are indicating that you want the state to be appended after the last existing state in the object.

## Important Commands To Know

### load

The "load" command is used to populate states of an object. By default, each new file loaded will be appended onto the object's states. However, the optional third argument to the load command is the frame index into which the file should be loaded. See "help load" or consult the reference section for more information.

```
load foo1.pdb,mov      # loads foo1.pdb into state 1 of "mov".
load foo2.pdb,mov      # loads foo2.pdb into state 2 of "mov".
load foo3.pdb,mov,3    # loads foo3.pdb into state 3 of "mov".
load foo4.pdb,mov,4    # loads foo4.pdb into state 4 of "mov".
```

### mset

The "mset" command is used to specify which states get included as frames of a movie. If the mset command is not used, PyMOL will by default play through all states in sequential order. However, if you wish to use the other movie commands (such as mdo), it is necessary to explicitly use the mset command to create a movie definition inside of PyMOL.

The mset command is followed by an arbitrarily list of statements which defines the entire movie. Each statement takes on one of three forms:

1. # A number simply indicates a state is to be played next.

2. **x#** A lowercase "x" immediately followed by a number (no space) indicates that the previous state should be repeated that many times total.
3. **–#** A hyphen immediately followed by a number (no space) indicates that a numeric sequence of states are to be appended onto the movie starting with the previously played state going to indicated state.

Once a movie has been defined, the red "VCR" controls in the lower–right–hand corner of the viewer can be used to step or play through the movie.

## Examples

```
mset 1 x30      # creates a 30 frame movie consisting of state 1 played 30 times.
mset 1 -30      # creates a 30 frame movie: states 1, 2, ... up through 30.
mset 1 -30 -2   # 58 frames: states 1, 2, ... , 29, 30, then 29, 28, ... , 4, 3, down to 2
mset 1 6 5 2 3  # 5 frames: states 1, 6, 5, 2, 3 in that order.
```

See "help mset" or the reference section for more information.

## mdo

The "mdo" command allows you to bind a particular series of PyMOL commands to a frame in the movie. For instance, you can perform a rotation about the axis at each frame of the movie in order to sweep the camera about the object. See "help mdo" or the reference section for more information.

**NOTE:** The "util" module includes two python commands for generating mdo commands, "util.mrock" and "util.mdo". These functions have not been documented, but the source code can be found in the file `modules/pymol/util.py`. Since they are actual python functions, explicit parenthesis are required to invoke them.

```
util.mrock(start, finish, angle, phase, loop-flag)
util.mroll(start, finish, loop-flag)
```

## mmatrix

The "mmatrix" command allows you to store and recall a particular viewing matrix to be used to set up frame 1 of the movie. This can be particularly helpful when you're trying to preserve a movie's orientation while performing other actions within PyMOL during the same session. See "help mmatrix" or the reference section for more information.

## Simple Examples

Here a static structure is subject to a gentle rock. The following statements create a sixty frame movie which simply rocks the protein by 10 degrees.

```
load test/dat/pept.pdb      # load a structure
mset 1 x60                  # define the movie
util.mrock(1,60,10,1,1)     # issues mdo commands to create +/- 10 deg. rock over 60 frames
```

In this next example, the protein is rotated through a full 360 sweep about the Y–axis over 120 frames

```
load test/dat/pept.pdb      # load a structure
mset 1 x120                  # define the movie
```

```
util.mroll(1,120,1)          # issues mdo commands to create full rotation over 120 frames
```

## Complex Examples

The following is a Python program (with a .py or .pym extension) which uses a Python loop to load a large number of numbered PDB files, and then configures PyMOL to show them both forwards and backwards.

```
from glob import glob
from pymol import cmd

file_list = glob("mov*.pdb"):

for file in file_list
    cmd.load(file,"mov")

cmd.mset("1 -%d -2"%len(file_list))
```

## Previewing Ray-traced Movie Images

PyMOL has the ability to cache a series of images in RAM and to play them back at a much higher rate than they could be rendered originally. This is most-useful for ray-traced images, but it can also be used with OpenGL images.

### cache\_frames

The **cache\_frames** option controls whether or not PyMOL saves frames in memory. Its usage is demonstrated in the following script. NOTE: caching images takes a tremendous amount of memory, so you should use the "viewport" command to shrink the window before utilizing this option.

```
viewport 320,240
load test/dat/pept.pdb
orient
hide
show sph
mset 1 x30
util.mrock 1,30,3,1,1
set ray_trace_frames=1
set cache_frames=1
mplay
```

### mclear

Once you have loaded a set of frames into RAM, the frames will remain there until you run the "mclear" command, even if you manipulate that model. You can also press the mclear button on the external GUI window.

```
mclear # flushes the frame cache
```

## Saving movies

## mpng

You can save movie images to numbered PNG format files with a common prefix. If you want each frame to be ray-traced, you should turn on raytracing of frames, turn off caching, and clear the cache (see the Movie Menu or use the following commands).

```
set ray_trace_frames=1
set cache_frames=0
mclear
```

You can save the movie using the "mpng" command, or you can save it from the "File" menu. Either way, you must provide a prefix which will be used to create numbered PNG files.

```
mpng mov # will create mov0001.png, mov0002.png, etc.
```

If you are compressing movies using Adobe Premiere (recommended for best quality), you will probably want to convert the files using ImageMagick or a similar package into a format that Premiere is capable of reading (such as ".tga" – targa format).



# Advanced Mouse Controls

## Picking Atoms and Bonds

The current mouse configuration is visible on the lower right-hand corner of the screen as a matrix. Under the default mouse configuration:

- In Viewing Mode:
  - ◆ A single atom can be "picked" using the "Pk1" function which is CTRL/right-click.
  - ◆ Multiple atoms (up to 4) can be "picked" in series using the "PkAt" function which is CTRL/right-middle.
- In Editing Mode:
  - ◆ A single bond can be "picked" using the "PkTB" function which is CTRL/right-click.
  - ◆ Multiple atoms (up to 4) can be "picked" in series using the "PkAt" function which is CTRL/right-middle.

Whenever an atom or bond has been picked, a number of atom selections are automatically defined as described in the following table:

- **(pk1)** The selected atom (or the first selected atom in a bond selection).
- **(pk2)** The second selected atom in a bond selection.
- **(pkfrag#)** A fragment of the molecule with its base adjacent to the selected bond or atom.
- **(pkchain)** The contiguous chain of atoms which contains the selected atom or bond.
- **(pkresi)** All atoms in the residue you picked.

You can click on a selection name to see visually which atoms are included in the selection. All of these selections can be used and manipulated as if they were created manually using the select command. Note however, that these selections are quite fragile, and will be automatically deleted in response to a number of common occurrences, such as loading a new object.

## Example Usage of the "pk" Atom Selections

Assuming that you picked an atom or bond...

```
show sticks,(pkresi)      # show sticks on the residue you picked

color read,(pkchain)      # color the chain you picked

remove (byres pk1)        # removes all atoms in the residue you picked
```

## The "lb" and "rb" Selections

Most of the time, the "pk1" atom selection will suffice. However, there are times when you need to specify two or more sets of atom selections. This is where "lb" and "rb" come in.

In addition to the "pk" atom selection set, there are two more atom selections which can be manipulated with the mouse. These are: (1) the left-button or "lb" selection and (2) the right-button or "rb" selection. Under the

default mouse configuration:

- The "(lb)" selection can be redefined using the "lb" function which is CTRL-SHIFT/left-click.
- The "(lb)" selection can be expanded using the "+lb" function which is CTRL/left-click.
- The "(rb)" selection can be redefined using the "rb" function which is CTRL-SHIFT/right-click.

Certain commands are designed to use "(lb)" and "(rb)" as default arguments. For instance, the "distance" command, if called without any arguments, will attempt to create a distance object between the (lb) and (rb) selections if they exist

```
# define (lb) by CTRL-SHIFT/left-clicking one atom
# define (rb) by CTRL-SHIFT/right-clicking another

dist # will create a distance object between the two atoms.
```

## Conformational Editing

Sorry, no documentation yet — these features won't be too useful until PyMOL is coupled up with an energy minimization engine.

# Crystallography Applications

## Crystal Symmetry

Ralf Grosse-Kunstleve has provided his SgLite module to enable PyMOL to deduce symmetry relationships from standard space group and unit cell information. Currently that information can only be provided to PyMOL as a CRYST1 record in the PDB file, which includes the correct space group identifier. However, it would be only a minor development task to add a means of assigning unit-cell and symmetry to any molecule object directly from the API.

The format of the CRYST1 record is as follows.

1 - 6	Record name	"CRYST1"	
7 - 15	Real(9.3)	a	a (Angstroms).
16 - 24	Real(9.3)	b	b (Angstroms).
25 - 33	Real(9.3)	c	c (Angstroms).
34 - 40	Real(7.2)	alpha	alpha (degrees).
41 - 47	Real(7.2)	beta	beta (degrees).
48 - 54	Real(7.2)	gamma	gamma (degrees).
56 - 66	LString	sGroup	Space group.
67 - 70	Integer	z	Z value. # ignored by PyMOL

## load

When you use the "load" command to read in a PDB file with symmetry information, matrix information should be output. Verify that this information is produced before attempting to display symmetry related molecules.

## symexp

The "symexp" command is used to display symmetry related molecules in the crystal lattice about an atom selection. This command creates a set of new objects with a common prefix. Each object in the series corresponds to one symmetry-related object, which can be treated independently. See "help symexp" or the reference section for usage information.

In order to visualize only symmetry-related atoms within a given distance, you need to break the process down into two steps. First, you use the symexp command to create complete symmetry-related objects. Then you use "hide" commands to restrict what is visible to only those areas which you are interested.

```
load foo.pdb                # load PDB file with CRYST record

symexp sym,foo,(foo),5.0    # Create symmetry related "foo" objects
                             # which pass within 5 angstroms of foo
                             # using the prefix "sym"

hide (not (foo expand 5))  # hide atoms greater than 5 A from foo
```

NOTE: The symexp command can potentially create large numbers of objects. You will want to use the "delete" command with a wildcard "\*" to remove all objects that share a common prefix.

```
delete sym*                # deletes objects starting with "sym"
```

# Electron Density Maps

The only map file format currently supported is the CNS and XPLOR ASCII format map file. PyMOL can read large maps of this format and then display arbitrary "bricks" of density within these maps about atom selections.

## load

PyMOL expects XPLOR/CNS map files to have a ".xplor" extension. This requirement can be avoided by supplying an explicit type of "xplor" to the "load" command.

```
load 2fofc.xplor,map1      # type inferred from the extension
load 2fofc.map,map1,1,xplor # type explicitly provided
```

See "help load" or the reference section for additional information.

## isomesh and isodot

Map objects are used to store the data and are represented by a wire-frame brick in space indicating the extent of the map. An arbitrary number of mesh or dots objects can be created from a given map using the "isomesh" and "isodot" commands.

```
isomesh msh1,map1,1.0 # display an isosurface-mesh at level 1.0 over
                      # the entire map object "map1"

isomesh msh2,map1,1.0,(chain A),3.0 # display isosurface-mesh at 1.0
                                     # in a brick about chain A with a
                                     # border of 3.0 Angstroms
```

See "help isomesh" or the reference section for additional information.

# Compiled Graphics Objects (CGOs) and Molscript Ribbons

## Introduction

Although PyMOL uses OpenGL for all real-time rendering, the simple ray-tracing engine inside of PyMOL is incapable of understanding arbitrary OpenGL calls. Thus, any graphics scene must be translated into a set of primitives (spheres, cylinders, and triangles) that can be provided to the ray-tracer in order to generate high quality images with "ideal" geometries, lighting, and shadows.

Compiled graphics objects are a PyMOL-specific format which enables any Python programmer to create three-dimensional geometries and animations which can be displayed at high-speed via OpenGL and also rendered into maximum-quality images via the raytracer without any additional work.

## Molscript Ribbons

NOTE: Molscript is a commercial software (free to academics) available at <http://www.avatar.se/molscript/> and must be obtained separately. It is our intention to eventually implement our own Molscript-quality ribbons directly from within PyMOL, but that day has not yet come.

PyMOL can automatically translate Raster3D format input files output by Molscript (with "-r" option) into Compiled Graphics Objects for display and rendering inside of PyMOL. PyMOL expects these files to have the file extension ".r3d". NOTE: the Raster3D-to-CGO interpreter is a bare-minimum Python implementation, and doesn't include anything beyond what is required to read what is output by Molscript.

## load

```
load test/dat/pept.r3d # loads one of the example raster3d files
```

## Using Molscript

### molauto

When using molauto to preparing input files for PyMOL, it is important to use the "-nocentre" option to prevent any transformation of the protein. That way the PDB file and the Molscript ribbons will be in the same frame of reference.

```
Unix> molauto -nocentre 3all.pdb | molscript -r > test1.r3d
Unix> molauto -nocentre -nice 3all.pdb | molscript -r > test2.r3d
```

You can load both PDB and ribbon files directly into PyMOL as separate objects.

```
load 3all.pdb # loads the coordinates
load test1.r3d # loads molscript ribbon
```

### Molscript Input Files

Unfortunately, PyMOL does not have the ability to write molscript input files which reflect the current atom

colorings and visibilities. Therefore, you will need to get in the habit of manually editing Molscript input files in order color and customize ribbons appropriately. Here are some tips:

1. Remove any line starting with "transform atom" from existing Molscript input files in order to preserve the frame of reference. For example:

```
transform atom * by centre position atom *;
```

2. For performance reasons, you may want to set the segments to a small number while working with Molscript ribbons in real-time. Later on you can increase this number, recreate, and reload the ".r3d" files.

```
set segments 2; # good for real-time graphics
```

```
set segments 8; # good for rendering
```

The easiest way to create new ribbons using PyMOL is to use the "save" command to write out a PDB file containing the atom selection of interest. You can then apply the "system" command to run molauto and molscript, and then load the Raster3D file back into PyMOL.

```
save tmp.pdb,(chain C)
system molauto -nocentre tmp.pdb | molscript -r > tmp.r3d
load tmp.r3d
```

## Creating Compiled Graphics Objects

Compiled graphics objects contain equivalents to the normal line and triangle primitives found in OpenGL but also include primitives for spheres and cylinders.

At the Python level, compiled graphics objects are constructed as a simple linear list of Python floating point numbers, which is conceptually equivalent to an OpenGL stream.

```
from pymol.cgo import *      # get constants
from pymol import cmd
```

```
obj = [
    BEGIN, LINES,
    COLOR, 1.0, 1.0, 1.0,

    VERTEX, 0.0, 0.0, 0.0,
    VERTEX, 1.0, 0.0, 0.0,

    VERTEX, 0.0, 0.0, 0.0,
    VERTEX, 0.0, 2.0, 0.0,

    VERTEX, 0.0, 0.0, 0.0,
    VERTEX, 0.0, 0.0, 3.0,

    END
]
```

```
cmd.load_cgo(obj, 'cgo01')
```

CGOs support the standard OpenGL BEGIN/END formalism as well as a few stand-alone primitives, SPHERE, CYLINDER, and TRIANGLE, which should NOT appear within a BEGIN/END block.

# CGO Reference

A CGO is simply a Python list of floating point numbers, which are compiled by PyMOL into a CGO object and associated with a given state.

Lowercase names below are should be replaced with floating-point numbers. Generally, the TRIANGLE primitive should only be used only as a last restore since it is much less effective to render than using a series of vertices with a BEGIN/END group.


```
BEGIN, { POINTS | LINES | LINE_LOOP | LINE_STRIP | TRIANGLES | TRIANGLE_STRIP | TRIANGLE_FAN },
VERTEX, x, y, z,
COLOR, red, green, blue,
NORMAL, normal-x, normal-y, normal-z,
END,
LINEWIDTH, line-width,
WIDTHSCALE, width-scale, # for ray-tracing
SPHERE, x, y, z, radius # uses the current color
CYLINDER, x1, y1, z1, x2, y2, z2, radius,
          red1, green1, blue1, red2, green2, blue2,
TRIANGLE, x1, y1, z1,
          x2, y2, z2,
          x3, y3, z3,
          normal-x1, normal-y1, normal-z1,
          normal-x2, normal-y2, normal-z2,
          normal-x3, normal-y3, normal-z3,
          red1, green1, blue1,
          red2, green2, blue2,
          red3, green3, blue3,
```

## load\_cgo

CGO lists are loaded into PyMOL using the "load\_cgo" function.

```
cmd.load_cgo(list,name,state)
```

Arbitrary 3D animations can be created by loading CGOs into consecutive states of a given object. The example below is a static image from the "examples/devel/cgo03.py" cgo animation demonstration program.

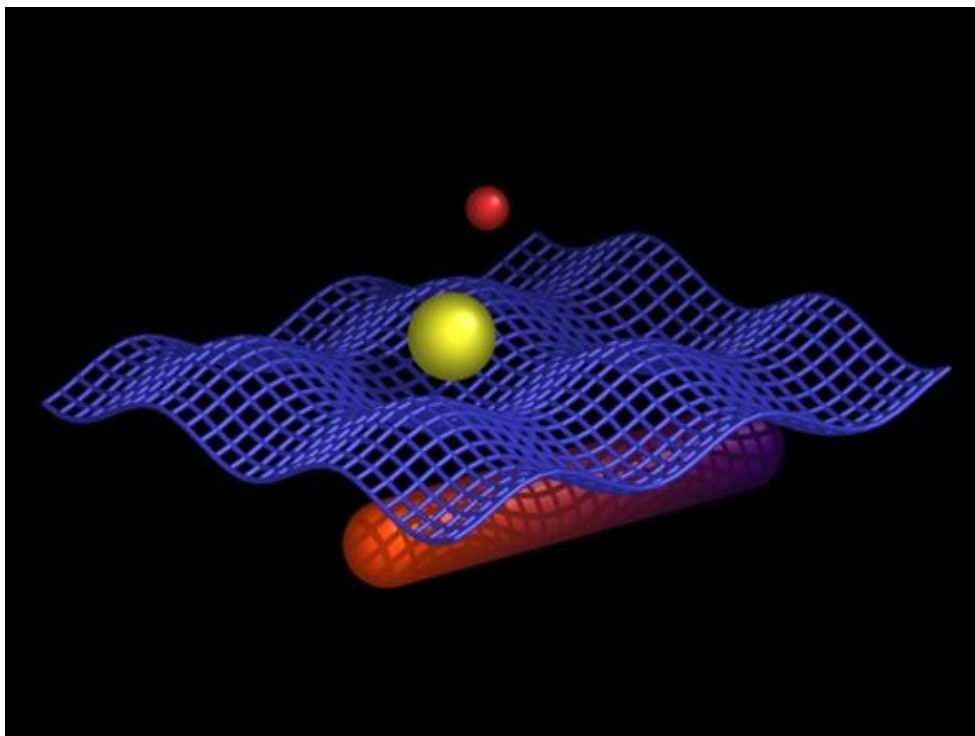



# Callback Objects and PyOpenGL

This is mainly a developer's function, so most users can skip this section. You will need some Python knowledge to understand the example.

## Introduction

Although all OpenGL rendering in PyMOL is performed at the C level, PyOpenGL provides an alternative binding of the OpenGL API from Python. Unfortunately, it is impossible for PyMOL to produce ray-traced images of objects rendered using PyOpenGL. Nevertheless, PyOpenGL can be used within PyMOL via Callback objects for pure OpenGL-based rendering purposes. If you need your graphics to be ray-traceable, then you should use Compiled Graphics Objects (see previous section).

Callback objects are trivial Python objects which have a "`__call__`" method for rendering and a "`get_extent`" method which tells PyMOL where in space the object is located. Once a callback object has been loaded into PyMOL, Python will automatically call this object when needed to update the screen.

PyMOL includes a copy of PyOpenGL under the pymol module hierarchy (`pymol.opengl`), but usage of this copy is of course optional. You can instead bind to the latest version without problems, provided that you install it yourself into the Python library that PyMOL is using by default.

NOTE: The current Windows version of PyMOL does not include Numeric, which makes heavy usage of PyOpenGL from within PyMOL impractical under Windows at present.

## Example

The following Python program shows how you can use a Callback object within PyMOL to perform rendering using OpenGL. For more examples, see the directory "`$PYMOL_PATH/examples/devcl`".

```
from pymol.opengl.gl import *
from pymol.callback import Callback
from pymol import cmd

class myCallback(Callback):

    def __call__(self):

        glBegin(GL_LINES)

        glColor3f(1.0,1.0,1.0)

        glVertex3f(0.0,0.0,0.0)
        glVertex3f(1.0,0.0,0.0)

        glVertex3f(0.0,0.0,0.0)
        glVertex3f(0.0,2.0,0.0)

        glVertex3f(0.0,0.0,0.0)
        glVertex3f(0.0,0.0,3.0)

        glEnd()

    def get_extent(self):
```

```
        return [[0.0,0.0,0.0],[1.0,2.0,3.0]]  
cmd.load_callback(myCallback(),'gl01')
```

## **load\_callback**

Callback objects are loaded into PyMOL using the "load\_callback" function.

```
cmd.load_callback(object,name,state)
```